



PRODUCT CATALOG

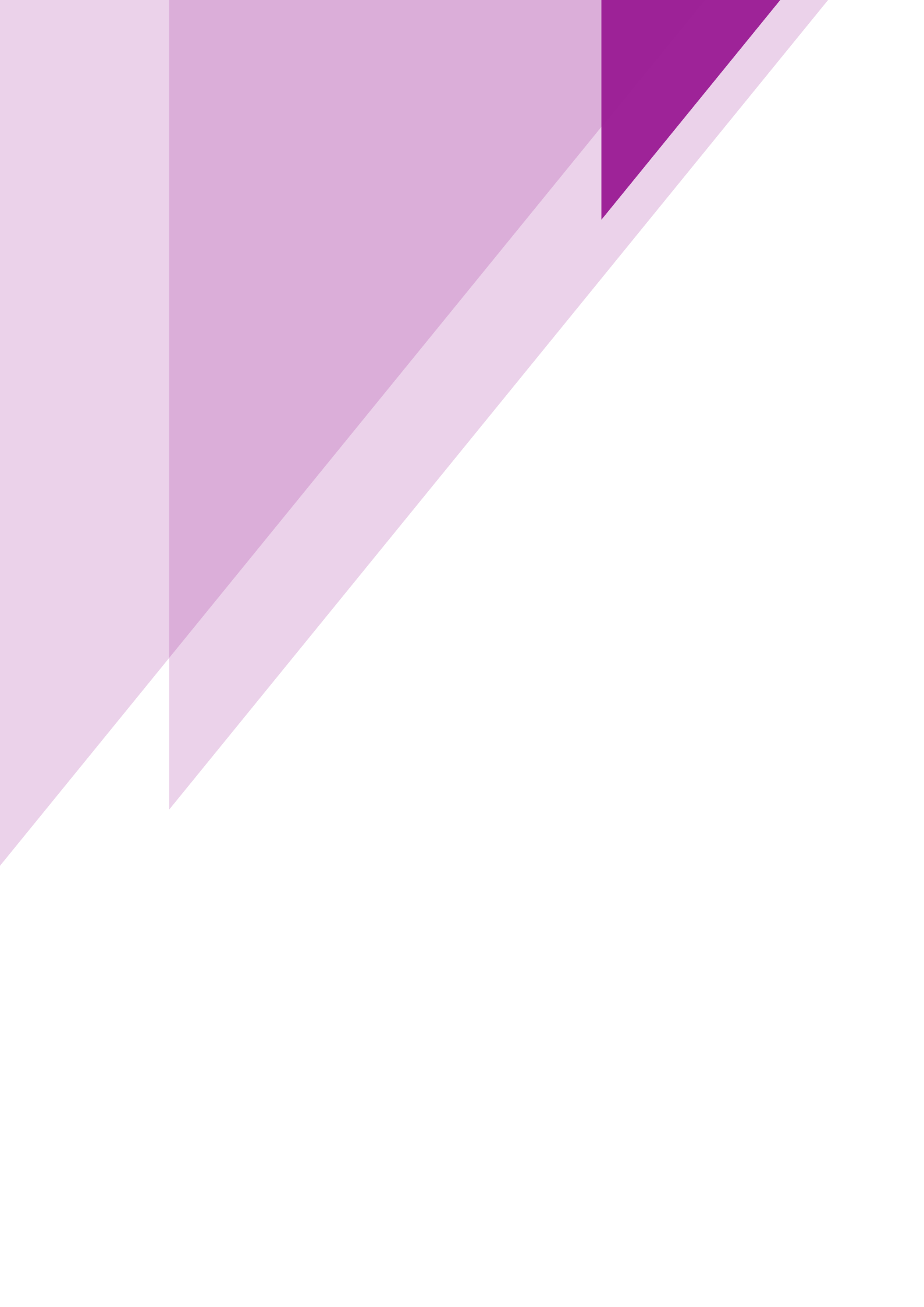



Table of contents

Standard Internal Hex Implant Systems	5	Fix-1™	46
Mars™	6 - 7	Ø 3.0	47
Ø 3.5	8	Ø 3.3	48
Ø 3.75	9	Ø 3.6	49
Ø 4.2	10	Ø 4.2	50
Ø 5.0	11	Ø 5.0	51
Ø 6.0	12	CAD CAM	52
FixTite™-S	13	Internal Hex	53-56
Ø 3.5	14	Fix Konisch™ ND	57 - 59
Ø 3.75	15	Fix Konisch™ PRO	60 - 62
Ø 4.2	16	Accessories	63 - 65
Ø 5.0	17	Standard Internal Hex Surgical Kit	66 - 67
Ø 6.0	18	Fix Konisch™ Surgical Kit	68 - 69
Fix-A-Dent™	19	Tools	70 - 71
Ø 3.3	20	Surgical Drills	72
Ø 3.75	21	Drill Stoppers	73
Ø 4.2	22	Guided Surgery Kit	74 - 77
Ø 5.0	23		
Ø 6.0	24		
Fix Konisch™ Conical Hex Connection	29 - 31		
Ø 3.0	32 - 36		
Ø 3.5	37		
Ø 4.3	38		
Ø 5.0	39 -45		

IVB

Standard Internal Hex Implant Systems


► Mars™

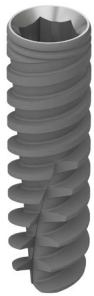

AB/AE - Alumina Oxide Blast / Acid Etched Surface Treatment



► FixTite™-S

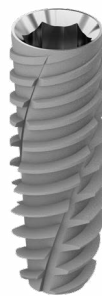

AB/AE - Alumina Oxide Blast/Acid Etched Surface Treatment



► Fix-a-Dent™



AB/AE - Alumina Oxide Blast / Acid Etched Surface Treatment



► Fix-1™



AB/AE - Alumina Oxide Blast / Acid Etched Surface Treatment




Surface Treatments

AB/AE (SLA™*) -

Alumina Oxide Blasted Large Grit Followed by Acid Etch

AB/AE surface treatment involves blasting with large grit aluminum oxide particles followed by acid etching with sulfuric and hydrochloric acid to remove the embedded particles.

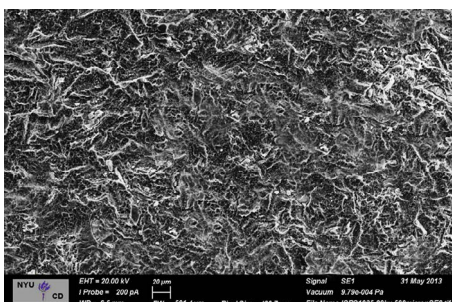

* Institut Straumann AG's trademark



Mars™


Mars™ dental implants are tapered implants with a spiral tap that condenses bone during placement for enhanced immediate stability. The Mars™ implants have two large variable threads and a tapered core for accurate implant placement, self-tapping, improved esthetics and better load distribution.

Mars™ implants feature a special round shaped apex that pushes the bone graft with minimal harm to anatomic structures. Mars™ implants are very user-friendly and their unique design often allows for immediate function for all bone types.

The Mars™ spiral dental implants feature a standard SLA-AB/AE- Alumina Oxide Blasted/Acid Etched surface treatment.

Mars™ features:

- High primary stability, also in type IV bone
- Double lead threads (2x1.2mm)
- Bone-condensing property
- Self-tapping
- Redirecting capability for optimal placement
- Built-in platform switching
- Minimally invasive surgery
- Immediate function / Loading capability
- Single prosthetic platform for easy restoration

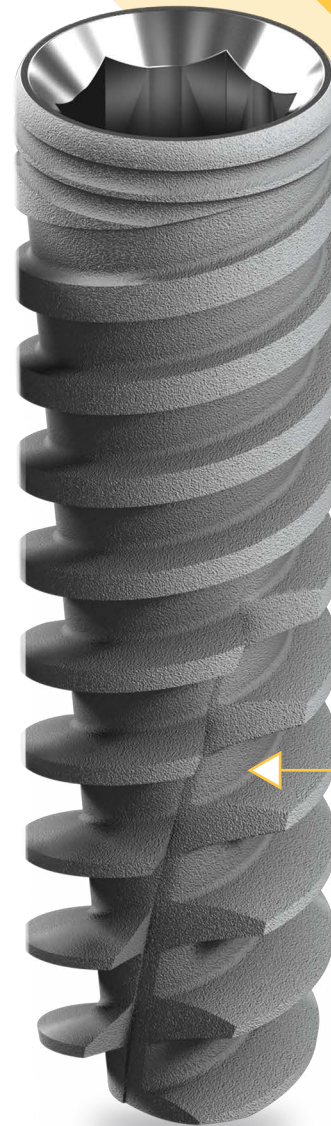


AB/AE Surface Treatment

	L-6.25mm	L-8mm	L-10mm	L-11.5mm	L-13mm	L-16mm	L-18mm
Ø3.5mm		113508	113510	113511	113513	113516	113518
Ø3.75mm		113808	113810	113811	113813	113816	113818
Ø4.2mm	114206	114208	114210	114211	114213	114216	114218
Ø5.0mm	115006	115008	115010	115011	115013	115016	
Ø6.0mm	116006	116008	116010	116011	116013		


Mars™ Ø 3.5

- ▶ AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

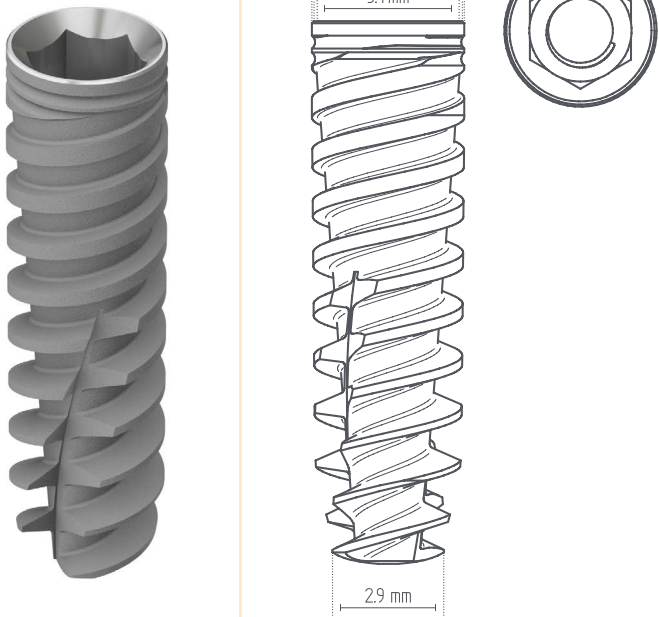Catalog No.	Length mm	
113508	8	
113510	10	
113511	11.5	
113513	13	
113516	16	
113518	18	

- ▶ Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step*	1. Tri-Step* 2. (3.2)	1. Tri-Step* 2. (3.2)

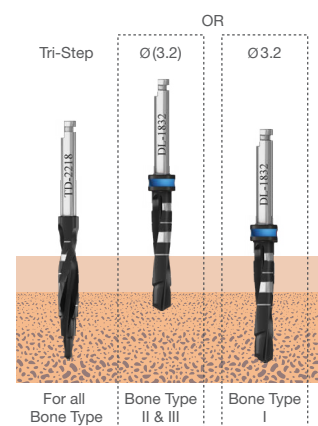



Note: (x.xx) drill to the depth of the cortex only. All measurements are in mm.

⚠ CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


Mars™ Ø 3.75

- ▶ AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

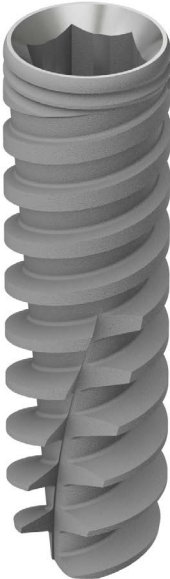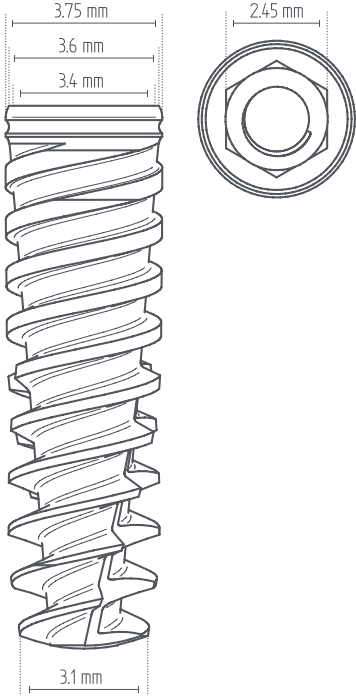Catalog No.	Length mm		
113808	8		
113810	10		
113811	11.5		
113813	13		
113816	16		
113818	18		

- ▶ Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step*	1. Tri-Step* 2. (3.2)	1. Tri-Step* 2. 3.2

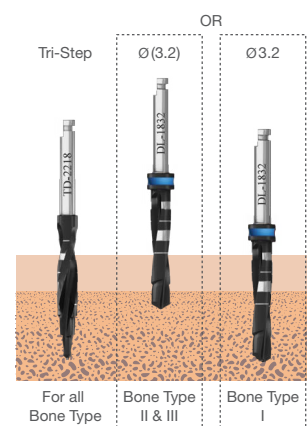


Note: (x.xx) drill to the depth of the cortex only. All measurements are in mm.

⚠ CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


Mars™ Ø 4.2

- ▶ AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

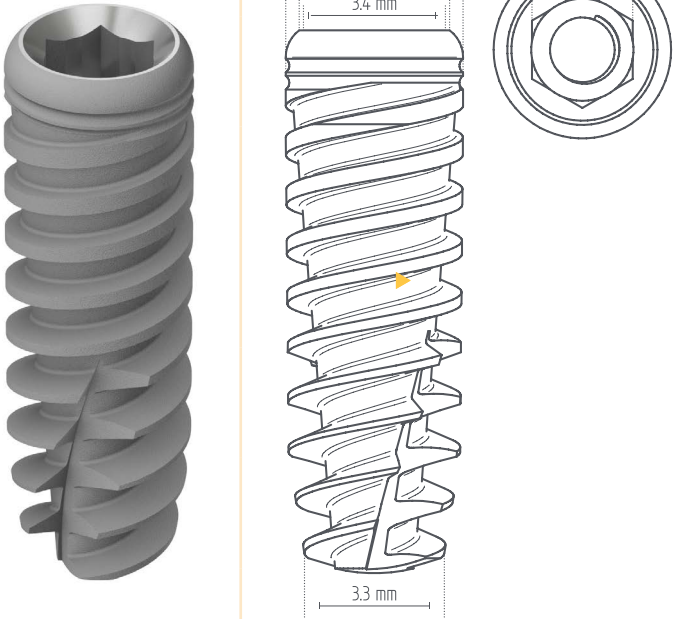Catalog No.	Length mm	
114206	6.25	
114208	8	
114210	10	
114211	11.5	
114213	13	
114216	16	
114218	18	

- ▶ Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step* 2. (3.6)	1. Tri-Step* 2. (3.6)	1. Tri-Step* 2. 3.6

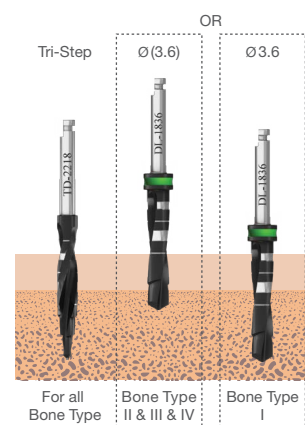


Note: (x.xx) drill to the depth of the cortex only. All measurements are in mm.

⚠ CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


Mars™ Ø 5.0

- ▶ AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

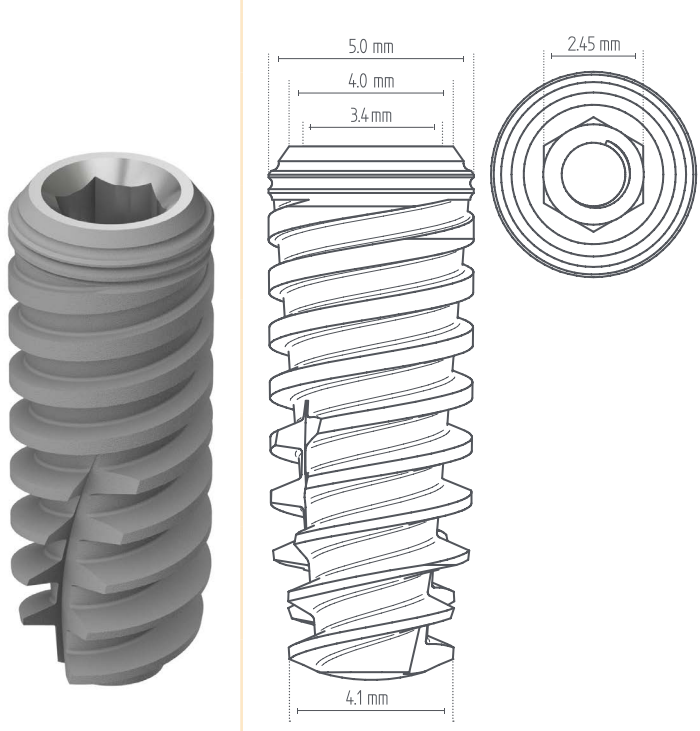Catalog No.	Length mm	
115006	6.25	
115008	8	
115010	10	
115011	11.5	
115013	13	
115016	16	

▶ Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step* 2. 3.6 3. (4.2)	1. Tri-Step* 2. 3.6 3. 4.2 4. (4.6)	1. Tri-Step* 2. 3.6 3. 4.2 4. 4.6

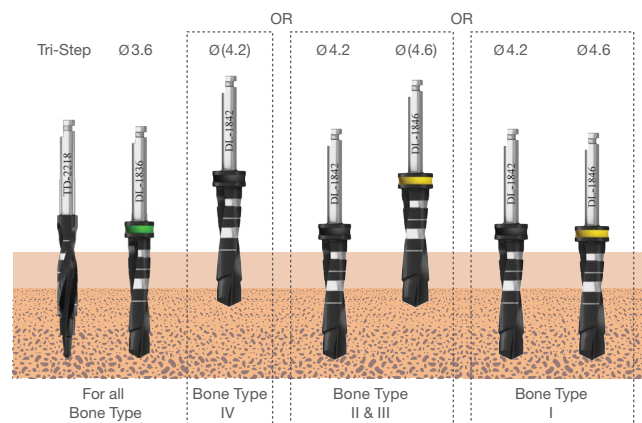

Note: (x.xx) drill to the depth of the cortex only.
All measurements are in mm.




CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


Mars™ Ø 6.0

- ▶ AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

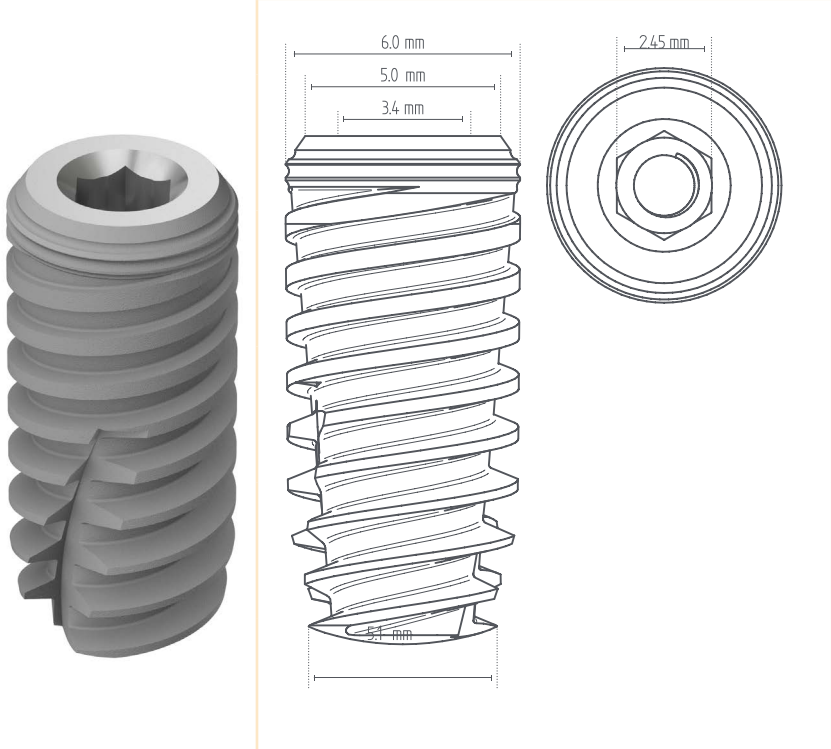Catalog No.	Length mm	
116006	6.25	
116008	8	
116010	10	
116011	11.5	
116013	13	

- ▶ Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step* 2. 3.6 3. 4.2 4. (5.2)	1. Tri-Step* 2. 3.6 3. 4.2 4. 5.2 5. (5.6)	1. Tri-Step* 2. 3.6 3. 4.2 4. 5.2 5. 5.6

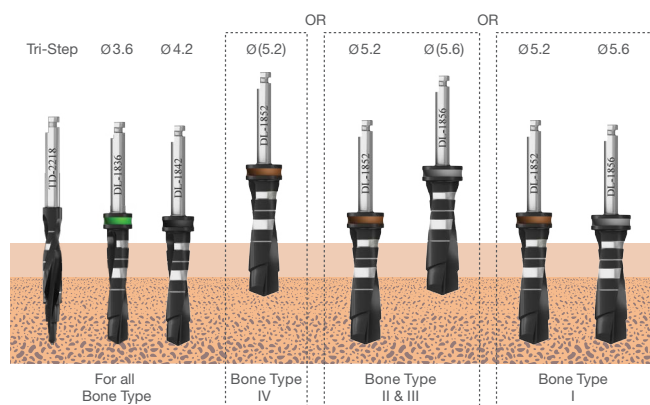

Note: (x.xx) drill to the depth of the cortex only.
All measurements are in mm.



CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.

FixTite™-S

FixTite™-S is a tapered core dental implant with a unique spiral triple-lead thread and triple vertical flutes design.

FixTite™-S has three major advantages:

Speed

FixTite™-S features a unique triple-lead thread design with 3.6mm (3x1.2mm) pitch. This design allows improved self-drilling, self-cutting and a very fast implant insertion.

Simplicity

A simple surgical procedure and drilling protocol based on three and two step drills sets the stage for accurate implant positioning. Each FixTite™-S implant features a 2.4mm internal hex connection regardless of its diameter. This single prosthetic platform significantly simplifies the restoration.


Stability

FixTite™-S has a conical core, spiral thread and triple vertical flutes design which creates immediate stability, also in type IV bone. The design contributes to balanced load distribution and reduces bone pressure with minimal invasive surgery.

FixTite™-S features a trapezoid thread and AB/AE surface treatment for accelerated osseointegration and increased BIC (Bone to Implant Contact).

FixTite™-S features:


- Triple-lead thread design with 3.6 (3x1.2mm) pitch
- Trapezoid thread, optimal space for a clot of blood
- Single restorative platform, regardless of diameter, with built-in platform switching
- AB/AE surface treatment
- Flat apex designed for improved insertion and placement
- Improved self-drilling and self-cutting
- Three vertical flutes

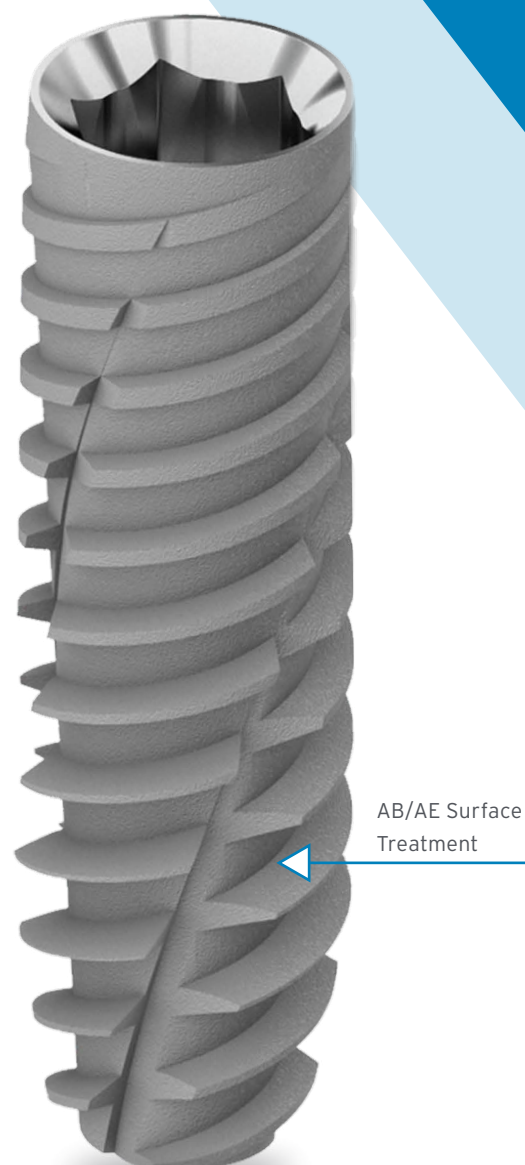

	L-6.25mm	L-8mm	L-10mm	L-11.5mm	L-13mm	L-16mm	L-18mm
Ø3.5mm		35081	35101	35111	35131	35161	35181
Ø3.75mm		37081	37101	37111	37131	37161	37181
Ø4.2mm	42061	42081	42101	42111	42131	42161	42181
Ø5.0mm	50061	50081	50101	50111	50131	50161	
Ø6.0mm	60061	60081	60101	60111	60131	60161	


FixTite™ - S Ø 3.5

- ▶ AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

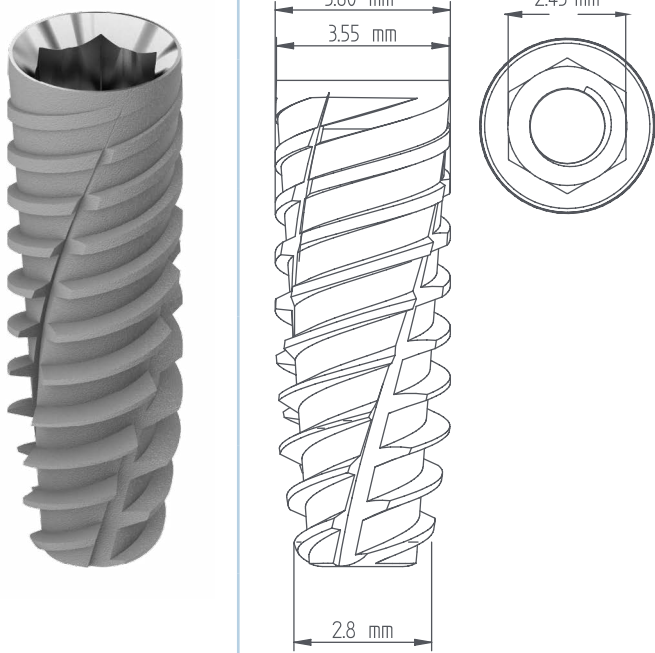Catalog No.	Length mm	
35081	8	
35101	10	
35111	11.5	
35131	13	
35161	16	
35181	18	

- ▶ Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step*	1. Tri-Step* 2. (3.2)	1. Tri-Step* 2. (3.2)

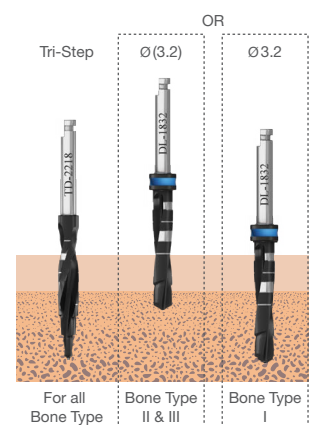

Note: (x.xx) drill to the depth of the cortex only. All measurements are in mm.




CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


FixTite™ - S Ø 3.75

- ▶ AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

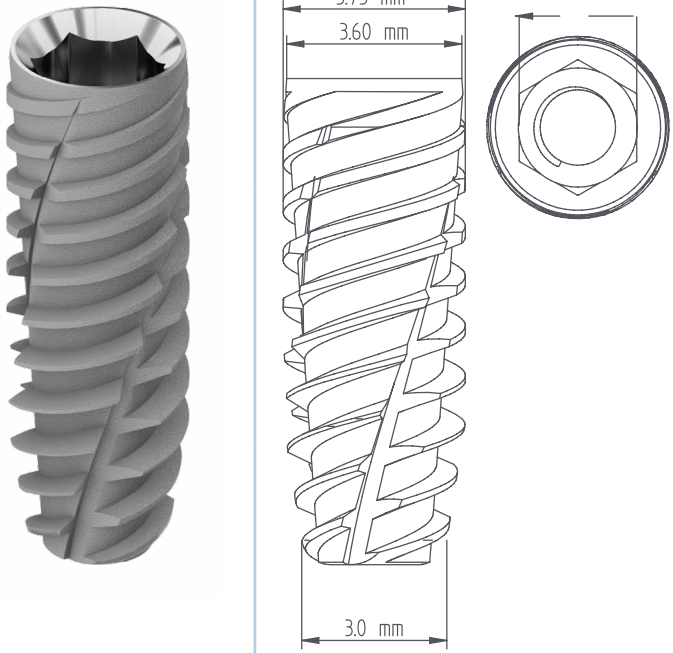Catalog No.	Length mm	
37081	8	
37101	10	
37111	11.5	
37131	13	
37161	16	
37181	18	

- ▶ Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step*	1. Tri-Step* 2. (3.2)	1. Tri-Step* 2. 3.2

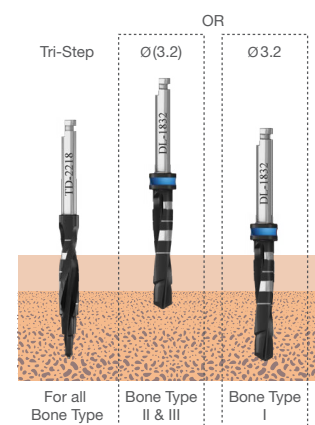


Note: (x.xx) drill to the depth of the cortex only. All measurements are in mm.

⚠ CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


FixTite™ - S Ø 4.2

- ▶ AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

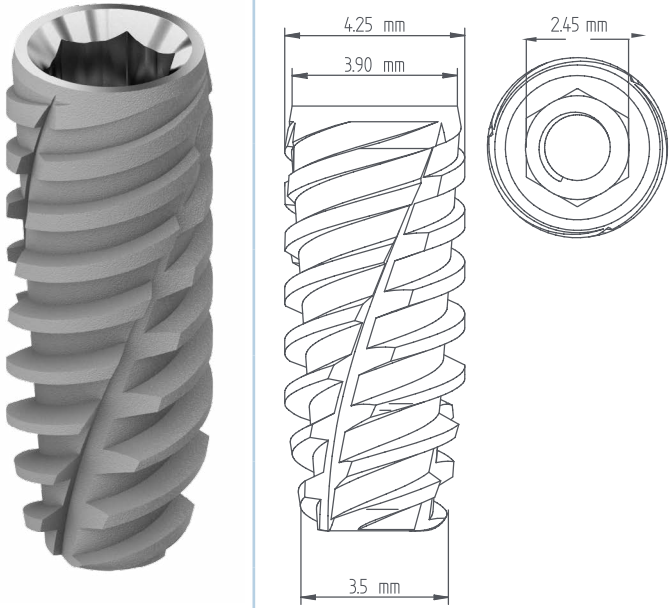Catalog No.	Length mm	
42061	6.25	
42081	8	
42101	10	
42111	11.5	
42131	13	
42161	16	
42181	18	

- ▶ Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step* 2. (3.6)	1. Tri-Step* 2. (3.6)	1. Tri-Step* 2. 3.6

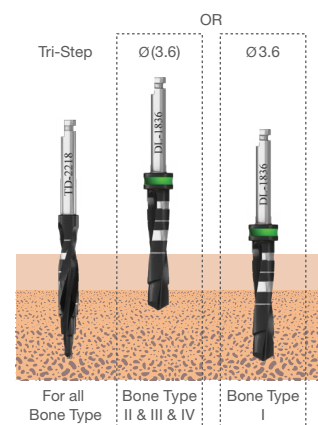

Note: (x.xx) drill to the depth of the cortex only. All measurements are in mm.




CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


FixTite™ - S Ø 5.0

- ▶ AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

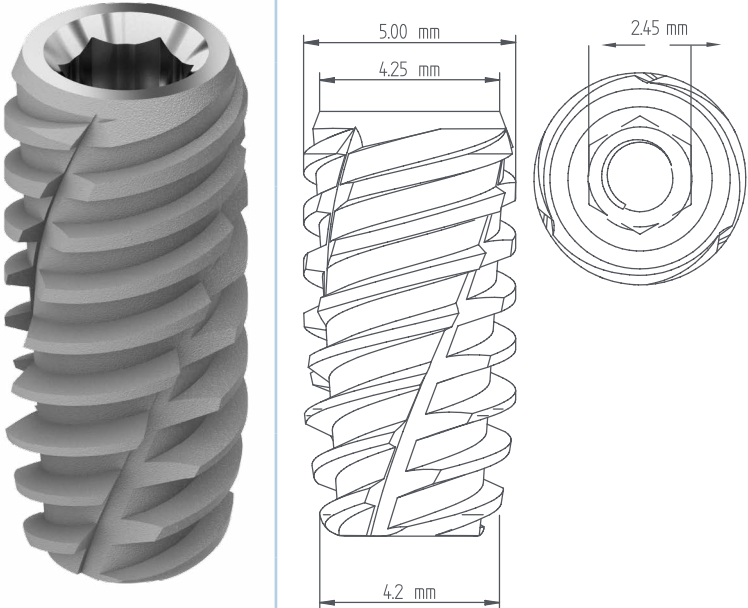Catalog No.	Length mm	
50061	6.25	
50081	8	
50101	10	
50111	11.5	
50131	13	
50161	16	

- ▶ Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step*	1. Tri-Step*	1. Tri-Step*
2. 3.6	2. 3.6	2. 3.6
3. (4.2)	3. 4.2	3. 4.2
	4. (4.6)	4. 4.6

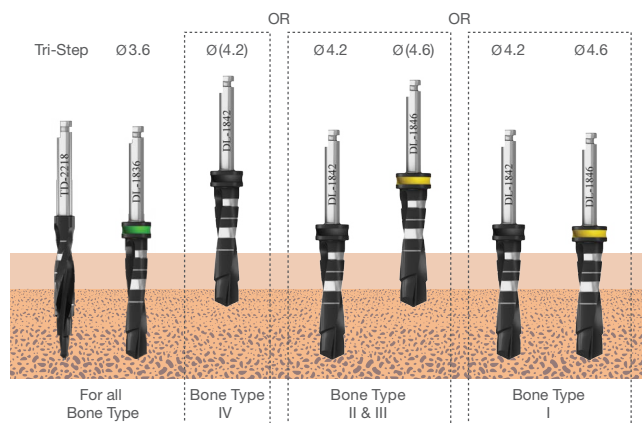

Note: (x.xx) drill to the depth of the cortex only.
All measurements are in mm.




CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


FixTite™ - S Ø 6.0

- ▶ AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

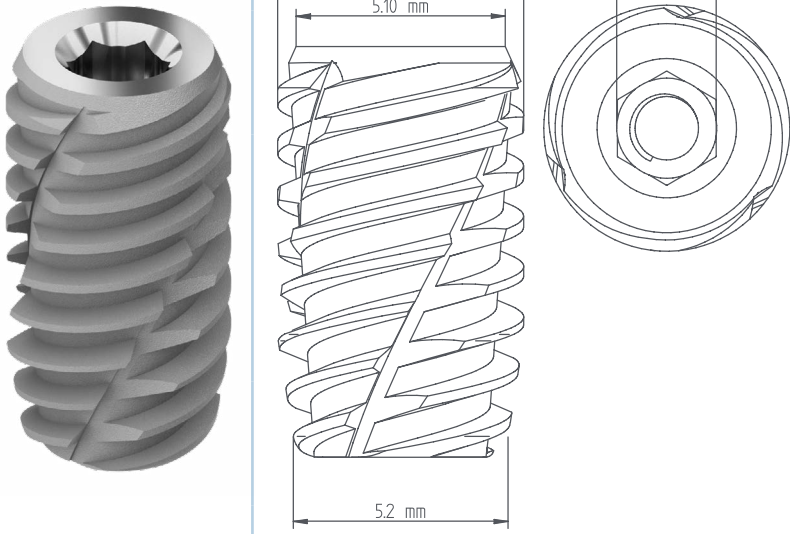Catalog No.	Length mm	
60061	6.25	
60081	8	
60101	10	
60111	11.5	
60131	13	
60161	16	

▶ Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step* 2. 3.6 3. 4.2 4. (5.2)	1. Tri-Step* 2. 3.6 3. 4.2 4. 5.2 5. (5.6)	1. Tri-Step* 2. 3.6 3. 4.2 4. 5.2 5. 5.6

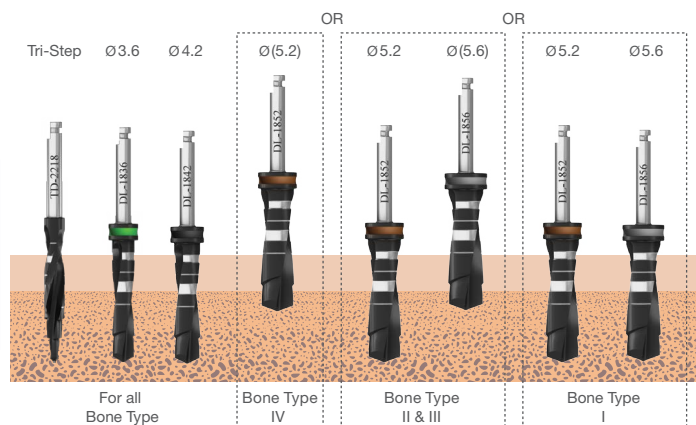

Note: (x.xx) drill to the depth of the cortex only.
All measurements are in mm.



CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.

Fix-A-Dent™

Fix-A-Dent™ Straight Implant System

The Fix-A-Dent™ implant is a straight, parallel-walled, slightly tapered implant, with a V-shaped thread design. Its design allows for accurate positioning to improve esthetics and to better load distribution. The Fix-A-Dent™ features AB/AE (Alumina-Oxide Blasted / Acid-Etched) surface treatment.


AB/AE Surface
Treatment



	L-6.25mm	L-8mm	L-10mm	L-11.5mm	L-13mm	L-16mm	L-18mm
Ø3.3mm			STA-3310	STA-3311	STA-3313	STA-3316	STA-3318
Ø3.75mm		STA-3708	STA-3710	STA-3711	STA-3713	STA-3716	STA-3718
Ø4.2mm	STA-4206	STA-4208	STA-4210	STA-4211	STA-4213	STA-4216	STA-4218
Ø5.0mm	STA-5006	STA-5008	STA-5010	STA-5011	STA-5013	STA-5016	
Ø6.0mm	STA-6006	STA-6008	STA-6010	STA-6011	STA-6013		


Fix-A-Dent™ - S Ø3.3

- AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

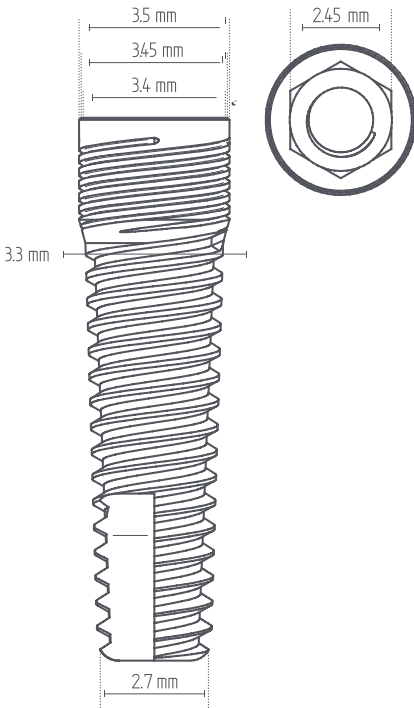Catalog No.	Length mm		
STA-3310	10		
STA-3311	11.5		
STA-3313	13		
STA-3316	16		
STA-3318	18		

- Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. 2.0 2. (2.8)	1. 2.0 2. (2.8) 3. (3.2)	1. Tri-Step*

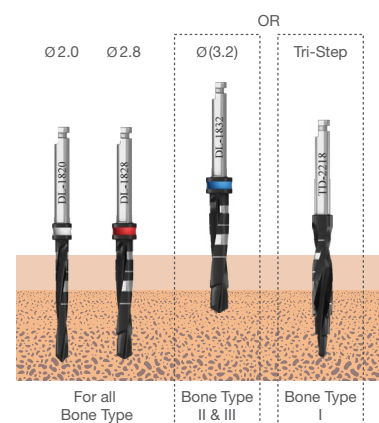

Note: (x.xx) drill to the depth of the cortex only.
All measurements are in mm.




CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


Fix-A-Dent™ - S Ø3.75

- AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

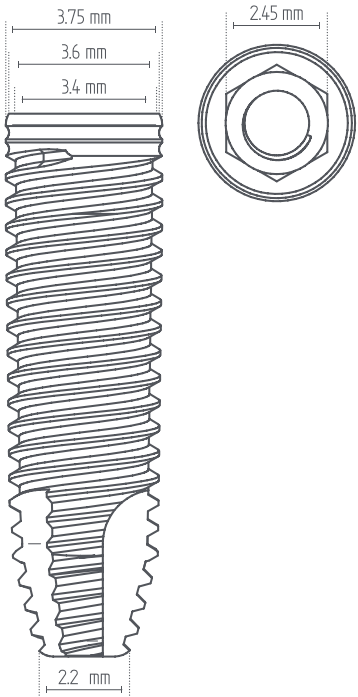Catalog No.	Length mm		
STA-3708	8		
STA-3710	10		
STA-3711	11.5		
STA-3713	13		
STA-3716	16		
STA-3718	18		

- Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step*	1. Tri-Step* 2. (3.2)	1. Tri-Step* 2. 3.2

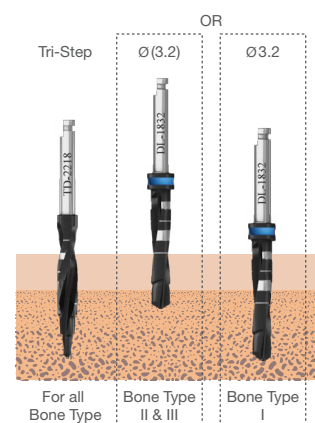


Note: (x.xx) drill to the depth of the cortex only. All measurements are in mm.

⚠ CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


Fix-A-Dent™ - S Ø4.2

- AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

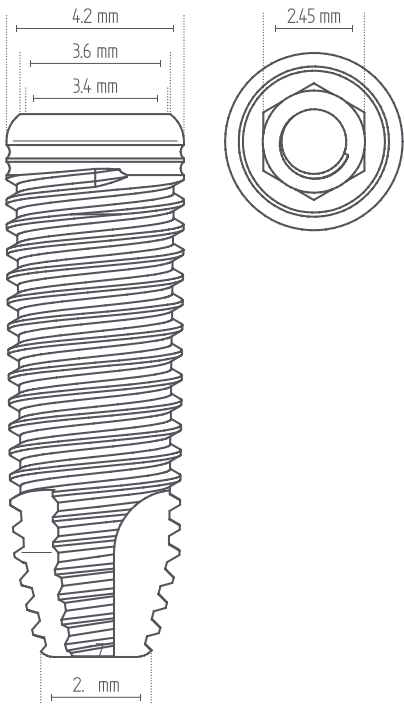Catalog No.	Length mm		
STA-4206	6.25		
STA-4208	8		
STA-4210	10		
STA-4211	11.5		
STA-4213	13		
STA-4216	16		
STA-4218	18		

- Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step* 2. 3.2	1. Tri-Step* 2. (3.6)	1. Tri-Step* 2. 3.6

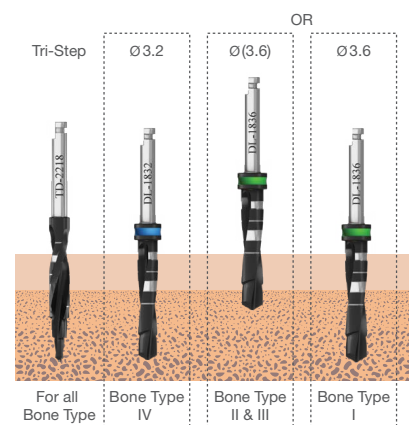

Note: (x.xx) drill to the depth of the cortex only.
All measurements are in mm.




CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


Fix-A-Dent™ - S Ø5.0

- AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

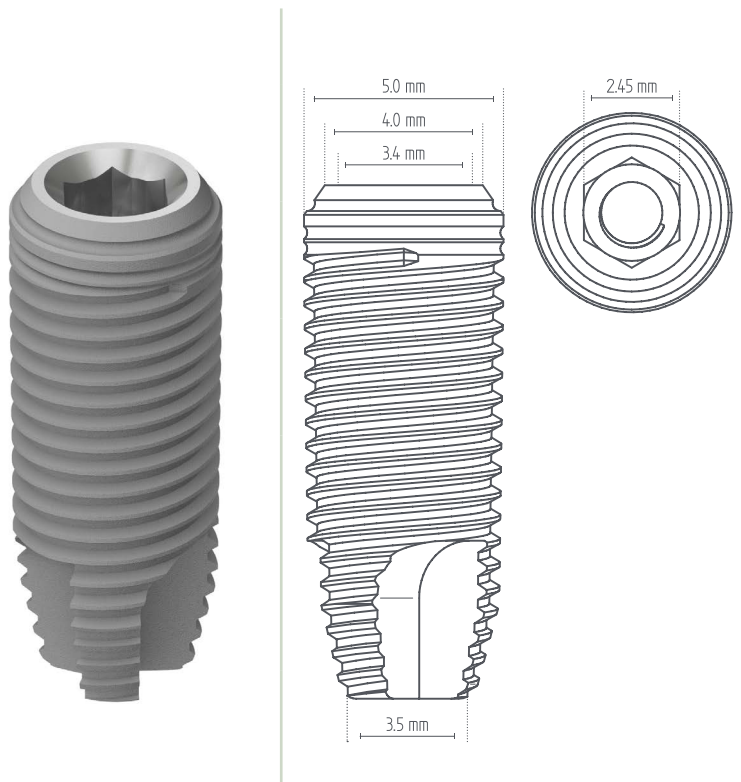Catalog No.	Length mm	
STA-5006	6.25	
STA-5008	8	
STA-5010	10	
STA-5011	11.5	
STA-5013	13	
STA-5016	16	

- Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step*	1. Tri-Step*	1. Tri-Step*
2. 3.6	2. 3.6	2. 3.6
3. (4.2)	3. 4.2	3. 4.2
	4. (4.6)	4. 4.6

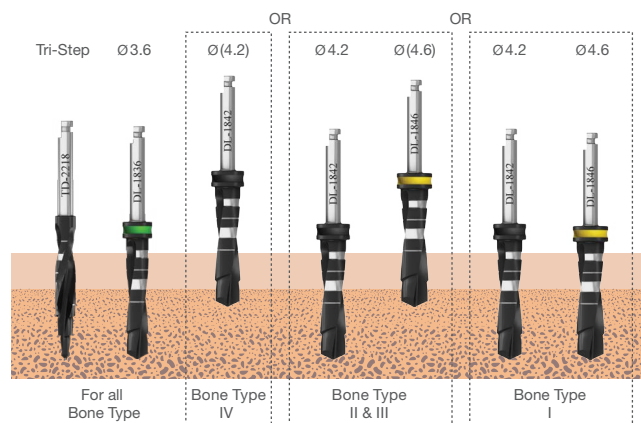


Note: (x.xx) drill to the depth of the cortex only.
All measurements are in mm.

⚠ CAUTION: The drill preparation is up to 1mm longer than the implant.
* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


Fix-A-Dent™ - S Ø6.0

- AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

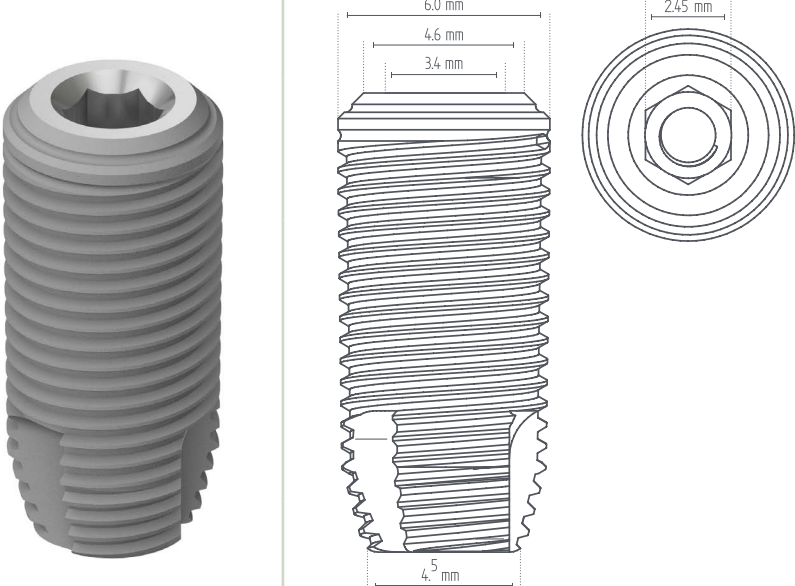Catalog No.	Length mm	
STA-6006	6.25	
STA-6008	8	
STA-6010	10	
STA-6011	11.5	
STA-6013	13	

- Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step*	1. Tri-Step*	1. Tri-Step*
2. 3.6	2. 3.6	2. 3.6
3. 4.2	3. 4.2	3. 4.2
4. 4.6	4. 4.6	4. 4.6
5. (5.2)	5. 5.2	5. 5.2
	6. (5.6)	6. (5.6)

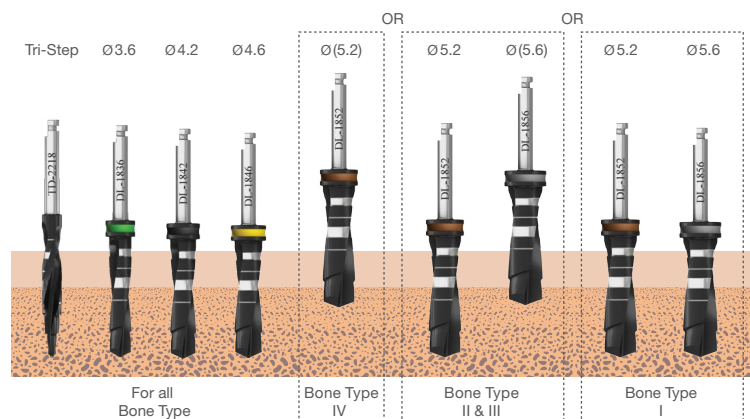

Note: (x.xx) drill to the depth of the cortex only.

All measurements are in mm.



CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.

Healing Abutments & Cover Screw

	CSS	Implant Cover Screw	Ti 6Al 4V ELI
	TGH2-SL	Slim Healing Abutment 3.5mmD x 2mmL	Ti 6Al 4V ELI
	TGH3-SL	Slim Healing Abutment 3.5mmD x 3mmL	Ti 6Al 4V ELI
	TGH4-SL	Slim Healing Abutment 3.5mmD x 4mmL	Ti 6Al 4V ELI
	TGH5-SL	Slim Healing Abutment 3.5mmD x 5mmL	Ti 6Al 4V ELI
	TGH6-SL	Slim Healing Abutment 3.5mmD x 6mmL	Ti 6Al 4V ELI
	TGH2-ST	Healing Abutment 4.5mmD x 2mmL	Ti 6Al 4V ELI
	TGH3-ST	Healing Abutment 4.5mmD x 3mmL	Ti 6Al 4V ELI
	TGH4-ST	Healing Abutment 4.5mmD x 4mmL	Ti 6Al 4V ELI
	TGH5-ST	Healing Abutment 4.5mmD x 5mmL	Ti 6Al 4V ELI
	TGH6-ST	Healing Abutment 4.5mmD x 6mmL	Ti 6Al 4V ELI
	TGH2-WD	Wide Healing Abutment 6mmD x 2mmL	Ti 6Al 4V ELI
	MIH3013	Wide Healing Abutment 6mmD x 3mmL	Ti 6Al 4V ELI
	TGH4-WD	Wide Healing Abutment 6mmD x 4mmL	Ti 6Al 4V ELI
	TGH5-WD	Wide Healing Abutment 6mmD x 5mmL	Ti 6Al 4V ELI

Cement Retained

	SRA3N-3	Slim Titanium Abutment 3mm	Ti 6Al 4V ELI
	SRA3N-6	Slim Titanium Abutment 6mm	Ti 6Al 4V ELI
	SRA3N-8	Slim Titanium Abutment 8mm	Ti 6Al 4V ELI
	SRA3	Straight Titanium Abutment - No End Line	Ti 6Al 4V ELI
	SRA3-1	Straight Titanium Abutment 1mm	Ti 6Al 4V ELI
	SRA3-2	Straight Titanium Abutment 2mm	Ti 6Al 4V ELI
	SRA3-3	Straight Titanium Abutment 3mm	Ti 6Al 4V ELI
	SRA3-4	Straight Titanium Abutment 4mm	Ti 6Al 4V ELI
	SRA3-5	Straight Titanium Abutment 5 mm	Ti 6Al 4V ELI
	TAL11.5	Straight Titanium Abutment Long 12.5mm	Ti 6Al 4V ELI
	SRA5	Straight Titanium Abutment 5mmD	Ti 6Al 4V ELI
	WPSRA	Wide Profile Titanium Abutment	Ti 6Al 4V ELI
	CAS2H	Straight Cylindrical Titanium Abutment 2mm - with Hex	Ti 6Al 4V ELI
	TAS2	Straight Cylindrical Titanium Abutment 2mm - without Hex	Ti 6Al 4V ELI

Cement Retained

	OPA	One Piece Screw-in Abutment	Ti 6Al 4V ELI
	PAB151	Angled Titanium Abutment 15°	Ti 6Al 4V ELI
	PAB151L	Angled Titanium Abutment 15° - Long	Ti 6Al 4V ELI
	PAB251	Angled Titanium Abutment 25°	Ti 6Al 4V ELI
	PAB351	Angled Titanium Abutment 35°	Ti 6Al 4V ELI
	PAB151-1	Esthetic 1mm Angled Titanium Abutment 15°	Ti 6Al 4V ELI
	PAB151-2	Esthetic 2mm Angled Titanium Abutment 15°	Ti 6Al 4V ELI
	PAB151-3	Esthetic 3mm Angled Titanium Abutment 15°	Ti 6Al 4V ELI
	CATPE	Titanium / Plastic Abutment with Hex	Ti 6Al 4V ELI/ POM
	CATPNE	Titanium / Plastic Abutment without Hex	Ti 6Al 4V ELI/ POM
	PRA3N	Plastic Cylinder Casting Abutment	POM
	PRA3N-NH	Plastic Cylinder Casting Abutment w/o Hex	POM
	PRA3	Plastic ST. Abutment Hex	POM
	PRA3-NH	Plastic Standard Abutment w/o Hex	POM
	PRA3R	Plastic Casting ST. Abutment Hex	POM
	PR5040	Plastic Cylinder Abutment Hex	POM

Screw Retained - MU Trans Mucosal Abutments

	TMI1	Straight Trans Mucosal Abutment 1mm	Ti 6Al 4V ELI
	TMI2	Straight Trans Mucosal Abutment 2mm	Ti 6Al 4V ELI
	TMI3	Straight Trans Mucosal Abutment 3mm	Ti 6Al 4V ELI
	TMI4	Straight Trans Mucosal Abutment 4mm	Ti 6Al 4V ELI
	TMI5	Straight Trans Mucosal Abutment 5mm	Ti 6Al 4V ELI
	TMA17	Angled Trans Mucosal Abutment 17° 2mm	Ti 6Al 4V ELI
	TMA1730	Angled Trans Mucosal Abutment 17° 3mm	Ti 6Al 4V ELI
	TMA30	Angled Trans Mucosal Abutment 30° 3mm	Ti 6Al 4V ELI
	TMA3040	Angled Trans Mucosal Abutment 30° 4mm	Ti 6Al 4V ELI
	TMA45	Angled Trans Mucosal Abutment 45° 4mm	Ti 6Al 4V ELI
	TMA-TGH	TMA Healing Cap	Ti 6Al 4V ELI

Screw Retained - MU Trans Muscosal Abutments

	TMA-TTC	TMA Titanium Temporary Cylinder	Ti 6Al 4V ELI
	TMA-PCS	TMA Plastic Casting Sleeve w / Prosthetic Screw	POM
	TMA-A	TMA Abutment Replica	Ti 6Al 4V ELI
	TMA-SRTO	TMA Impression Copying Open Tray	Ti 6Al 4V ELI
	TMA-SRTC	TMA Closed Tray Transfer	Ti 6Al 4V ELI

Ball Attachment - Overdentures

	TMIB1-P3	Ball Attachment 1mm w / Hex	Ti 6Al 4V ELI
	TMIB2-P3	Ball Attachment 2mm w / Hex	Ti 6Al 4V ELI
	TMIB3-P3	Ball Attachment 3mm w / hex	Ti 6Al 4V ELI
	TMIB4-P3	Ball Attachment 4mm w / Hex	Ti 6Al 4V ELI
	TMIB5-P3	Ball Attachment 5mm w / Hex	Ti 6Al 4V ELI
	TMIB6-P3	Ball Attachment 6mm w / Hex	Ti 6Al 4V ELI
	NC01	Plastic Ball Cap - White	Nylon
	NCP	Plastic Ball Cap Pink	Nylon
	HOT	Ball Cap Titanium	Ti 6Al 4V ELI

Analog

	FIA	Internal Hex Implant Analog	Ti 6Al 4V ELI
---	-----	-----------------------------	---------------

Transfer / Impression Coping

	SOT10	Open Tray Slim Transfer 3.6mmDx10mmL	Ti 6Al 4V ELI
	SOT14	Open Tray Slim Transfer 3.6mmDx14mmL	Ti 6Al 4V ELI
	SCT09	Slim Closed Tray Transfer 3.5mmDx9mmL	Ti 6Al 4V ELI
	SCT13	Slim Closed Tray Transfer 3.5mmDx13mmL	Ti 6Al 4V ELI

* Transfers will be chosen by length and width according to the planned prosthetic treatment and according to the interocclusal distance.

IH Retrieval Screw

	SRAS1	Retrieval Screw for Abutments	Ti 6Al 4V ELI
---	-------	-------------------------------	---------------

IH Screws

	CSS	Implant Cover Screw	Ti 6Al 4V ELI
	SRAS	Prosthetic Screw	Ti 6Al 4V ELI
	SRAS2	Short Prosthetic Screw	Ti 6Al 4V ELI
	PAS	Screw for Angled Titanium Abutment	Ti 6Al 4V ELI

IH Screw			
	SRAS3	Retaining Screw For Angled TMA	Ti 6Al 4V ELI
	TPS	TMA Prosthetic Screw	Ti 6Al 4V ELI
	TMA-SRAS1	Retaining Screw for TMA 1mm	Ti 6Al 4V ELI
	TMA-SRAS2	Retaining Screw for TMA 2mm	Ti 6Al 4V ELI
	TMA-SRAS3	Retaining Screw for TMA 3mm	Ti 6Al 4V ELI
	TMA-SRAS4	Retaining Screw for TMA 4mm	Ti 6Al 4V ELI
	TMA-SRAS5	Retaining Screw for TMA 5mm	Ti 6Al 4V ELI
	TMA-SRAS-OT	TMA Screw for Open Tray Impression Copying	Ti 6Al 4V ELI
	CTSS	Closed Tray Transfer Screw Short	Ti 6Al 4V ELI
	CTSL	Closed Tray Transfer Screw Long	Ti 6Al 4V ELI

* All systems are provided with their compatible screws. However, they may also be ordered individually.

Prosthetic Hand Drivers			
	SHHS	1.27 Hex Hand Driver - Short	Stainless Steel
	SHHM	1.27 Hex Hand Driver - Medium	Stainless Steel
	SHHL	1.27 Hex Hand Driver - Long	Stainless Steel
	TMA-SHHL	TMA/FC Female Hand Hex Driver - Long	Stainless Steel
	TMA-SHHS	TMA/FC Female Hand Hex Driver - Short	Stainless Steel

Torque Drivers			
	SMTS	2.4 Hex Torque Driver - Short	Stainless Steel
	SMTL	2.4 Hex Torque Driver - Long	Stainless Steel
	SHT	1.27 Hex Torque Driver - Short	Stainless Steel
	SHTL	1.27 Hex Torque Driver - Long	Stainless Steel
	TDBA	Torque Driver for Ball Attachment	Stainless Steel
	TMA-SHT	TMA / FC Hex Torque Driver - Short	Stainless Steel
	TMA-SHTL	TMA / FC Hex Torque Driver - Long	Stainless Steel

Handpiece Drivers and Adaptors			
	TMA-HAS	TMA / FC Hex Handpiece Adapter - Short	Stainless Steel
	TMA-HAL	TMA / FC Hex Handpiece Adapter - Long	Stainless Steel
	SHHA9021	RS Hex Handpiece Adapter Short	Stainless Steel
	MHHA9026	RS Hex Handpiece Adapter Medium	Stainless Steel
	LHHA9029	RS Hex Hand Piece Adapter Long	Stainless Steel

Fix Konisch™

Fix Konisch™ PRO

Ø3.5mm

Standard Diameter (SD)
for Regular Size Ridges

Fix Konisch™ PRO

5.0mm

Wide Diameter (WD)
for Wide / Large Ridges



ND Fix Konisch™

Ø3.0mm

Narrow Diameter (ND)
Ridges and Tight Spaces

Fix Konisch™ PRO

Ø4.3 Wide Diameter (WD)
for Wide / Large Ridges

Fix Konisch™ Conical Hex Connection

The Fix Konisch™ spiral implant system offers a unique, strong and solid conical-hex connection. The conical-hex connection enables implants to be used in all indications and regions in the lower and upper jaw. The Morse Taper Connection minimizes micro movements between the implant and abutment.

The Fix Konisch™ spiral implant system is comprised of a tapered core implant with a spiral tap, promoting increased immediate stability. With their unique tip and double lead thread design (2x1.2mm), these implants slice through bone unlike conventional self-tapping implants, which scrape the bone away as they tap.

By repositioning bone as it is being inserted, Fix Konisch™ implants can be placed at different angles for exceptional primary stability and optimal restorative orientation, along with improved esthetics and a better loading distribution.

All Fix Konisch™ implants feature the RBM surface treatment proven to achieve the desired roughness levels needed for optimal osseointegration, increase the rate of Bone to Implant Contact (BIC) and attain the highest implant surface purity levels.

Fix Konisch™ Advantages:

- RBM Surface Treatment
- High Primary Stability, Even in Type IV Bone
- Self-Tapping
- Double Lead Threads (2 X 1.2mm)
- Redirecting Capability for more Accurate Implant Placement
- Built-In Platform Switching
- Strong Fix Konisch™ Connection
- Requires Minimally Invasive Surgery
- Immediate Function / Loading Capability in all Bone Types

	L-8mm	L-10mm	L-11.5mm	L-13mm	L-15mm	L-16mm	L-18mm
ND-Ø3.0mm		FKND1030	FKND1130	FKND1330		FKND1630	FKND1830
PRO SD-Ø3.5mm	FKSW0835	FKSW1035	FKSW1135	FKSW1335	FKSW1535		FKSW1835
PRO WD-Ø4.3mm	FKSW0843	FKSW1043	FKSW1143	FKSW1343	FKSW1543		FKSW1843
PRO WD-Ø5.0mm	FKSW0849	FKSW1049	FKSW1149	FKSW1349	FKSW1549		FKSW1849

Fix Konisch™ Implant System

Conical Hex Connection with RBM

Calcium Phosphate Blasted Surface Treatment



Actual color of implant platforms

Packaging implant diameter color code

RBM

RBM treatment process is achieved by “high-speed particle blasting”, using suitable & resorbable hardened Calcium Phosphate particles. Usage of Calcium Phosphate as the blasting material also eliminates the need of using strong acids for the removal of blasting material remnants.

Following the mechanical abrasion, leftover calcium particles are cleaned off by etching the surface with organic acid of low concentration, leaving no residue and avoiding implant’s titanium surface refurbishment.



RBM Mechanical & Biological Properties:

Successful osteogenesis around the implant, is obtained due to created homogenous pore diameter and uniform surface. Calcium Phosphate is a resorbable, biocompatible substance that does not transfix deeply into the implants surface and does not lodge inside the physically created crater.

In case there are Calcium Phosphate leftover residue on the implant surface, it adds to the osteogenic value of the implant due to the osteoconductive properties it possesses, allowing it to be completely resorbed and replaced by new bone cells. Implant surfaces that have been RBM treated are considered to be more osteoconductive. Research showed higher percentage of bone to implant contact around RBM treated implants. This finding could be especially useful in severe clinical conditions like poor bone quality and in cases of early or immediate loading.


RBM Surface Treatment Advantages:

- Cleaner Implant surface
- Mineralized bone formation contributed by optimal homogenous pore diameter.
- No Titanium nor acid residues left.
- Accelerated Healing & Osseointegration.
- Enlarged implant surface area for optimal contact

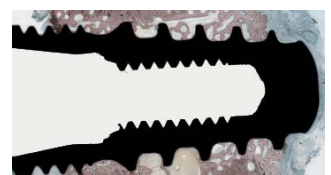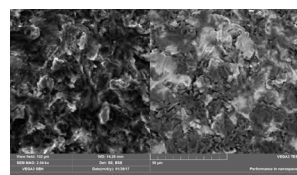


Fix-Konisch™ Narrow (ND) Dental Implant Ø3.0

- RBM Surface
Titanium Alloy - Ti 6Al 4V ELI

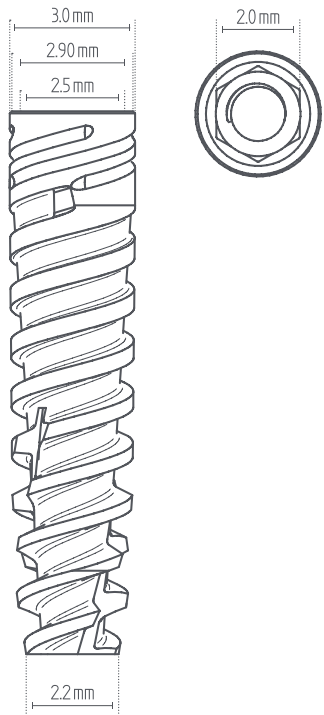Catalog No.	Length mm		
FKND1030	10		
FKND1130	11.5		
FKND1330	13		
FKND1630	16		
FKND1830	18		

► Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. 2.0	1. 2.0 2. (2.8)	1. 2.0 2. 2.8

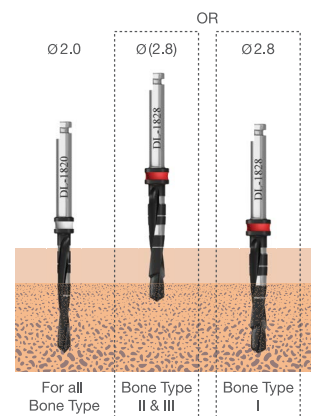

Note: (x.xx) drill to the depth of the cortex only. All measurements are in mm.

⚠ CAUTION: The drill preparation is up to 1mm longer than the implant.

Healing Abutments & Cover Screw

	MND0073	ND Conical Implant Cover Screw	Ti 6Al 4V ELI
	MND0046	ND Slim Healing Abutment 3mmD x 2mmL	Ti 6Al 4V ELI
	MND0047	ND Slim Healing Abutment 3mmD x 3mmL	Ti 6Al 4V ELI
	MND0048	ND Slim Healing Abutment 3mmD x 4mmL	Ti 6Al 4V ELI
	MND0049	ND Slim Healing Abutment 3mmD x 5mmL	Ti 6Al 4V ELI
	MND0050	ND Slim Healing Abutment 3mmD x 6mmL	Ti 6Al 4V ELI
	MND0004	ND Healing Abutment 4.5mmD x 2mmL	Ti 6Al 4V ELI
	MND0008	ND Healing Abutment 4.5mmD x 3mmL	Ti 6Al 4V ELI
	MND0009	ND Healing Abutment 4.5mmD x 4mmL	Ti 6Al 4V ELI
	MND0012	ND Healing Abutment 4.5mmD x 5mmL	Ti 6Al 4V ELI
	MND0013	ND Healing Abutment 4.5mmD x 6mmL	Ti 6Al 4V ELI
	MND0090	ND Conical Healing Abutment 4.5mmD x 2mmL	Ti 6Al 4V ELI
	MND0091	ND Conical Healing Abutment 4.5mmD x 3mmL	Ti 6Al 4V ELI
	MND0092	ND Conical Healing Abutment 4.5mmD x 4mmL	Ti 6Al 4V ELI
	MND0093	ND Conical Healing Abutment 4.5mmD x 5mmL	Ti 6Al 4V ELI
	MND0094	ND Conical Healing Abutment 4.5mmD x 6mmL	Ti 6Al 4V ELI

Cement Retained

	MND0010	ND Straight Slim Abutment - No End Line	Ti 6Al 4V ELI
	MND0019	ND Straight Abutment 1mm	Ti 6Al 4V ELI
	MND0011	ND Straight Abutment 2mm	Ti 6Al 4V ELI
	MND0064	ND Straight Abutment 3mm	Ti 6Al 4V ELI
	MND0018	ND Angled Abutment 15° - No End Line	Ti 6Al 4V ELI
	MND0015	ND Angled Abutment 15° 1mm	Ti 6Al 4V ELI
	MND0026	ND Angled Abutment 15° 2mm	Ti 6Al 4V ELI
	MND0063	ND Angled Abutment 15° 3mm	Ti 6Al 4V ELI
	MND0025	ND Angled Abutment 25° - No End Line	Ti 6Al 4V ELI
	MND0043	ND Anti-rotation Temporary Abutment	Ti 6Al 4V ELI
	MND0043	ND Free Rotation Temporary Abutment	Ti 6Al 4V ELI
	MND0044	ND Plastic Cylinder Abutment	POM

Long Laboratory Screw Flat Connection

	MND3095	ND Anti-rotation Titanium/Plastic Abutment	Ti 6Al 4V ELI / POM
	MND3097	ND Free-rotation Titanium/Plastic Abutment	Ti 6Al 4V ELI / POM

Screw Retained - Multi Unit Abutment (MU)

	MND0032	ND Straight Multi Unit Abutment 1mm	Ti 6Al 4V ELI
	MND0033	ND Straight Multi Unit Abutment 2mm	Ti 6Al 4V ELI
	MND0034	ND Straight Multi Unit Abutment 3mm	Ti 6Al 4V ELI
	MND3729	ND Straight Multi Unit Abutment 4mm	Ti 6Al 4V ELI
	MND3730	ND Straight Multi Unit Abutment 5mm	Ti 6Al 4V ELI
	MND0035	ND Angled Multi Unit Abutment °17 L=3mm	Ti 6Al 4V ELI
	MND0038	ND Angled Multi Unit Abutment °17 L=4mm	Ti 6Al 4V ELI
	MND0036	ND Angled Multi Unit Abutment °30 L=3.5mm	Ti 6Al 4V ELI
	MND0037	ND Angled Multi Unit Abutment °30 L=5mm	Ti 6Al 4V ELI
	TMA-TGH	TMA Healing Cap	Ti 6Al 4V ELI
	TMA-TTC	TMA Titanium Temporary Copyingr	Ti 6Al 4V ELI
	TMA-PCS	TMA Plastic Casting Sleeve w/Prosthetic Screw	POM
	TMA-A	TMA Abutment Replica	Stainless Steel
	TMA-SRTO	TMA Impression Copying Open Tray	Ti 6Al 4V ELI
	TMA-SRTC	TMA Closed Tray Transfer	Ti 6Al 4V ELI

* The term "impression coping" has been amended to the term "transfer" (the terms are synonymous). Please note that either term may appear in different documentation.

Screw Retained - Flat Connection

	MND0030	ND Flat Connection Abutment 2mm	Ti 6Al 4V ELI
	MND0031	ND Flat Connection Abutment 4mm	Ti 6Al 4V ELI
	FLAT5515	Temporary Cylinder Sleeve Flat Connection	Ti 6Al 4V ELI
	FLAT5737	Implant Replica Flat Connection	Stainless Steel
	FLAT6011	Castable Abutment for Centering Ring Flat Connection	POM
	FLAT6015	Titanium Retaining Screw Flat Connection	Ti 6Al 4V ELI
	FLAT6013	Titanium Gluing Ring Flat Connection	Ti 6Al 4V ELI
	FLAT6017	Long Laboratory Screw Flat Connection	Ti 6Al 4V ELI
	FLAT6018	Centering Ring Flat Connection	Ti 6Al 4V ELI
	FLAT5512	Transfer Flat Connection	Ti 6Al 4V ELI

Bell Attachment - Overdentures			
	MND0020	ND Ball Attachment 0.5mm	Ti 6Al 4V ELI
	MND0021	ND Ball Attachment 1mm	Ti 6Al 4V ELI
	MND0022	ND Ball Attachment 2mm	Ti 6Al 4V ELI
	MND0023	ND Ball Attachment 3mm	Ti 6Al 4V ELI
	MND0024	ND Ball Attachment 4mm	Ti 6Al 4V ELI
	NC01	Plastic Ball Cap white	Nylon
	NCP	Plastic Ball Cap Pink Soft	Nylon
	HOT	Ball Cap Titanium	Ti 6Al 4V ELI
Analog			
	MND0007	ND Implant Analog	Ti 6Al 4V ELI
	MIH5740	Ball Attachment Abutment Replica	Stainless Steel
Transfers/Impression Coping			
	MND5113	ND Fix Konisch™ Slim Close Tray Transfer 3.2mmDx13mmL	Ti 6Al 4V ELI
	MND5213	ND Fix Konisch™ Close Tray Transfer 4.2mmDx13mmL	Ti 6Al 4V ELI
	MND5614	ND Fix Konisch™ Slim Open Tray Transfer 3.2mmDx14mmL	Ti 6Al 4V ELI
	MND5514	ND Fix Konisch™ Open Tray Transfer 4.2mmDx14mmL	Ti 6Al 4V ELI
ND Screws			
	MND0073	ND Conical Implant Cover Screw	Ti 6Al 4V ELI
	MND0001	ND Abutment Screw	Ti 6Al 4V ELI
	MND0039	ND Retaining Screw for MU 1mm	Ti 6Al 4V ELI
	MND0040	ND Retaining Screw for MU 2mm	Stainless Steel
	MND0041	ND Retaining Screw for MU 3mm	POM
	MND3409	ND Retaining Screw for MU 4mm	Ti 6Al 4V ELI
	MND3412	ND Retaining Screw for MU 5mm	Ti 6Al 4V ELI
	MND0042	ND Retaining Screw for Angled MU	Ti 6Al 4V ELI
	MND5014	ND Fix Konisch™ Open Tray Transfer Screw	Ti 6Al 4V ELI
	MND5127	ND Fix Konisch™ Close Tray Transfer Screw	Ti 6Al 4V ELI
	MND0014	ND Retrieval Screw for Abutments	Ti 6Al 4V ELI

* All systems are provided with their compatible screws. However, they may also be ordered individually.

Hand Drivers

	MSD0064	Star Hand Driver-Short	Stainless Steel
	MSD0004	Star Hand Driver - Long	Stainless Steel
	TMA-SHHS	TMA/FC Female Hand Hex Driver - Short	Stainless Steel
	TMA-SHHL	TMA/FC Female Hand Hex Driver - Long	Stainless Steel

Torque Drivers


	MND1001	ND Abutment Insertion Tool	Stainless Steel
	MND0006	ND 2mm Hex Torque Driver - Short	Stainless Steel
	MND0005	ND 2mm Hex Torque Driver - Long	Stainless Steel
	MSD0065	Star Torque Driver-Short	Stainless Steel
	MSD0005	Star Torque Driver - Long	Stainless Steel
	TDBA	Torque Driver for Ball Attachment	Stainless Steel
	TMA-SHT	TMA/FC Hex Torque Driver - Short	Stainless Steel
	TMA-SHTL	TMA/FC Hex Torque Driver - Long	Stainless Steel
	WRS	Wrench Ratchet Surgical Instruments Handle	Stainless Steel

Handpiece Drivers and Adaptors

	MND0068	ND 2mm Hex Handpiece Adapter - Short	Ti 6Al 4V ELI
	MND0071	ND 2mm Hex Handpiece Adapter - Medium	Ti 6Al 4V ELI
	MND0074	ND 2mm Hex Handpiece Adapter - Long	Ti 6Al 4V ELI
	MSD0018	Star Handpiece Driver - Short	Ti 6Al 4V ELI
	MSD0019	ND 2mm Hex Handpiece Adapter - Short	Stainless Steel
	TMA-HAS	ND 2mm Hex Handpiece Adapter - Medium	Stainless Steel
	TMA-HAL	ND 2mm Hex Handpiece Adapter - Long	Stainless Steel


Fix-Konisch™ PRO Standard Dental Implant Ø3.5

- RBM Surface
Titanium Alloy - Ti 6Al 4V ELI

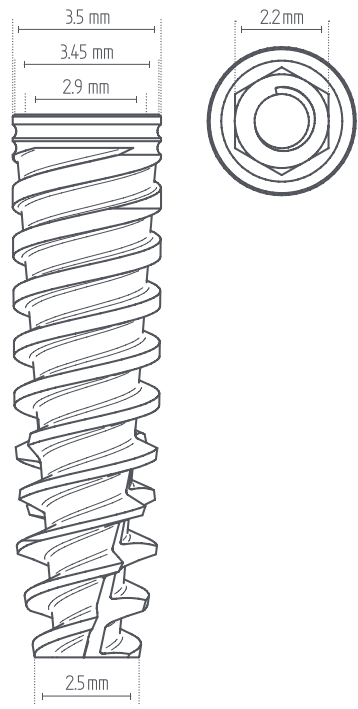Catalog No.	Length mm		
FKSW0835	8		
FKSW1035	10		
FKSW1135	11.5		
FKSW1335	13		
FKSW1535	15		
FKSW1835	18		

- Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step*	1. Tri-Step* 2. (3.2)	1. Tri-Step* 2. (3.2)

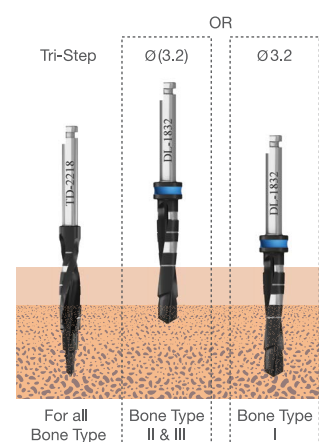



Note: (x.xx) drill to the depth of the cortex only. All measurements are in mm.

⚠ CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


Fix-Konisch™ PRO Wide Dental Implant Ø4.3

- RBM Surface
Titanium Alloy - Ti 6Al 4V ELI

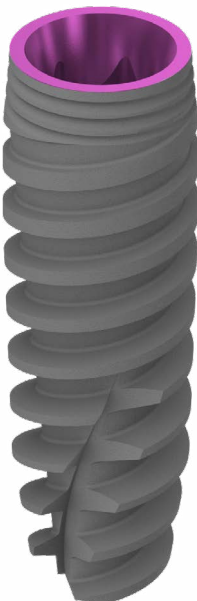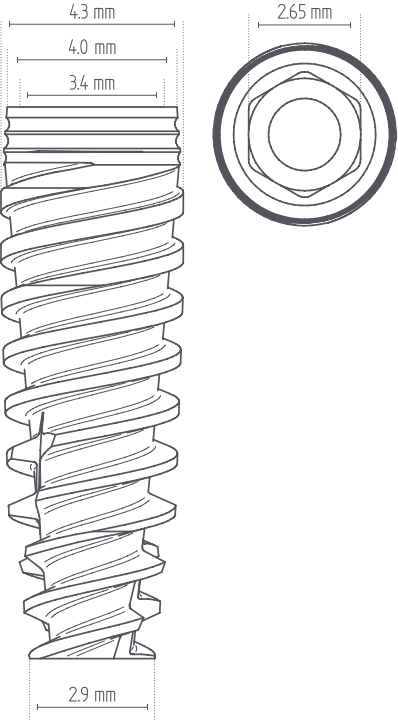Catalog No.	Length mm		
FKSW0843	8		
FKSW1043	10		
FKSW1143	11.5		
FKSW1343	13		
FKSW1543	15		
FKSW1843	18		

- Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step* 2. (3.6)	1. Tri-Step* 2. (3.6)	1. Tri-Step* 2. 3.6

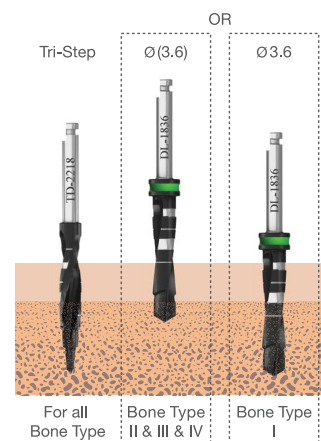



Note: (x.xx) drill to the depth of the cortex only. All measurements are in mm.

⚠ CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


Fix-Konisch™ PRO Wide Dental Implant Ø5.0

- RBM Surface
Titanium Alloy - Ti 6Al 4V ELI

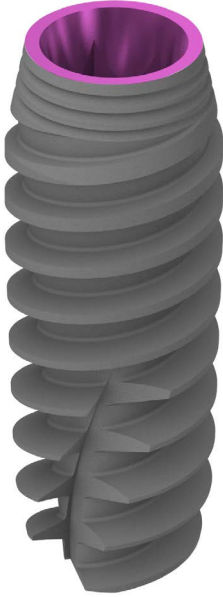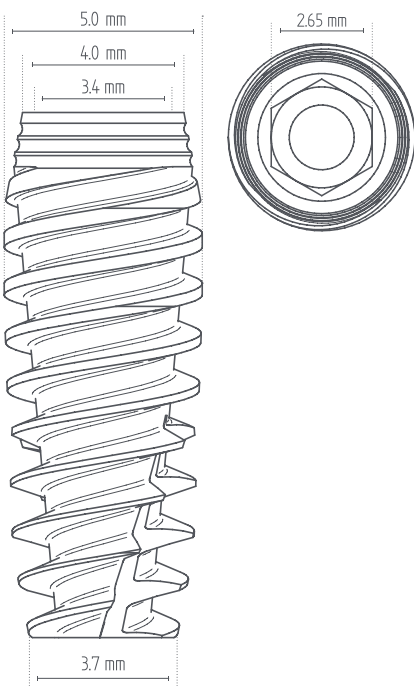Catalog No.	Length mm		
FKSW0849	8		
FKSW1049	10		
FKSW1149	11.5		
FKSW1349	13		
FKSW1549	15		
FKSW1849	18		

- Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step*	1. Tri-Step*	1. Tri-Step*
2. 3.6	2. 3.6	2. 3.6
3. (4.2)	3. 4.2	3. 4.2
	4. (4.6)	4. 4.6

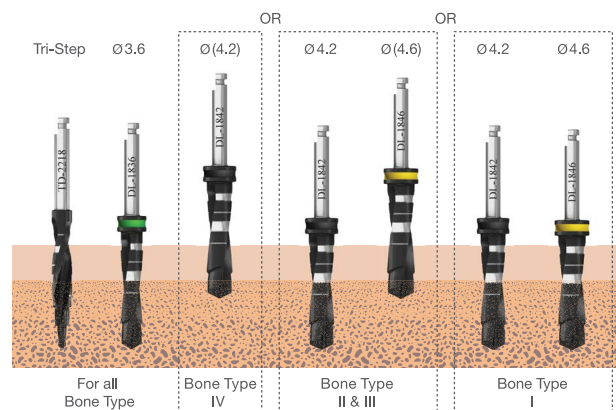

Note: (x.xx) drill to the depth of the cortex only.
All measurements are in mm.



CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.

Healing Abutments & Cover Screw

	MSD0073	SD Conical Implant Cover Screw	Ti 6Al 4V ELI
	MSD0046	SD Slim Healing Abutment 3.5mmD x 2mmL	Ti 6Al 4V ELI
	MSD0047	SD Slim Healing Abutment 3.5mmD x 3mmL	Ti 6Al 4V ELI
	MSD0048	SD Slim Healing Abutment 3.5mmD x 4mmL	Ti 6Al 4V ELI
	MSD0049	SD Slim Healing Abutment 3.5mmD x 5mmL	Ti 6Al 4V ELI
	MSD0050	SD Slim Healing Abutment 3.5mmD x 6mmL	Ti 6Al 4V ELI
	MSD0022	SD Healing Abutment 4.5mmD x 2mmL	Ti 6Al 4V ELI
	MSD0023	SD Healing Abutment 4.5mmD x 3mmL	Ti 6Al 4V ELI
	MSD0024	SD Healing Abutment 4.5mmD x 4mmL	Ti 6Al 4V ELI
	MSD0025	SD Healing Abutment 4.5mmD x 5mmL	Ti 6Al 4V ELI
	MSD0026	SD Healing Abutment 4.5mmD x 6mmL	Ti 6Al 4V ELI
	MSD0090	SD Conical Healing Abutment 4.5mmD x 2mmL	Ti 6Al 4V ELI
	MSD0091	SD Conical Healing Abutment 4.5mmD x 3mmL	Ti 6Al 4V ELI
	MSD0092	SD Conical Healing Abutment 4.5mmD x 4mmL	Ti 6Al 4V ELI
	MSD0093	SD Conical Healing Abutment 4.5mmD x 5mmL	Ti 6Al 4V ELI
	MSD0094	SD Conical Healing Abutment 4.5mmD x 6mmL	Ti 6Al 4V ELI

Cement Retained

	MSD0020	SD Straight Slim Abutment - No End Line	Ti 6Al 4V ELI
	MSD0012	SD Straight Abutment 1mm	Ti 6Al 4V ELI
	MSD0011	SD Straight Abutment 2mm	Ti 6Al 4V ELI
	MSD0010	SD Straight Abutment 3mm	Ti 6Al 4V ELI
	MSD0013	SD Angled Abutment °15 - No End Line	Ti 6Al 4V ELI
	MSD0014	SD Angled Abutment 1°15mm	Ti 6Al 4V ELI
	MSD0015	SD Angled Abutment 2°15mm	Ti 6Al 4V ELI
	MSD0043	SD Angled Abutment 3°15mm	Ti 6Al 4V ELI
	MSD0044	SD Angled Abutment °25 - No End Line	Ti 6Al 4V ELI
	MSD0045	SD Anti-rotation Temporary Abutment	POM
	MSD3095	SD Free Rotation Temporary Abutment	Ti 6Al 4V ELI
	MSD3097	SD Plastic Cylinder Abutment	Ti 6Al 4V ELI

Screw Retained - Multi Unit Abutment (MU)

	MSD0032	SD Straight Multi Unit Abutment 1mm	Ti 6Al 4V ELI
	MSD0033	SD Straight Multi Unit Abutment 2mm	Ti 6Al 4V ELI
	MSD0034	SD Straight Multi Unit Abutment 3mm	Ti 6Al 4V ELI
	MSD3729	SD Straight Multi Unit Abutment 4 mm	Ti 6Al 4V ELI
	MSD3730	SD Straight Multi Unit Abutment 5 mm	Ti 6Al 4V ELI
	MSD3734	SD Angled Multi Unit Abutment 2.5 °17mm	Ti 6Al 4V ELI
	MSD3737	SD Angled Multi Unit Abutment 3.5 °17mm	Ti 6Al 4V ELI
	MSD3735	SD Angled Multi Unit Abutment 3 °30mm	Ti 6Al 4V ELI
	MSD3738	SD Angled Multi Unit Abutment 4 °30mm	Ti 6Al 4V ELI
	TMA-TGH	TMA Healing Cap	Ti 6Al 4V ELI
	TMA-TTC	TMA Titanium Temporary Copying	Ti 6Al 4V ELI
	TMA-PCS	TMA Plastic Casting Sleeve w/Prosthetic Screw	POM
	TMA-A	TMA Abutment Replica	Ti 6Al 4V ELI
	TMA-SRTO	TMA Impression Copying Open Tray	Ti 6Al 4V ELI
	TMA-SRTO	MA Closed Tray Transfer	Ti 6Al 4V ELI

* The term "impression coping" has been amended to the term "transfer" (the terms are synonymous). Please note that either term may appear in different documentation.

Screw Retained - Flat Connection

	MSD0030	SD Abutment Flat Connection 2mm	Ti 6Al 4V ELI
	MSD0031	SD Abutment Flat Connection 4mm	Ti 6Al 4V ELI
	FLAT5515	Temporary Cylinder Sleeve Flat Connection	Ti 6Al 4V ELI
	FLAT5737	Implant Replica Flat Connection	Stainless Steel
	FLAT6011	Castable Abutment for Centering Ring Flat Connection	POM
	FLAT6015	Titanium Retaining Screw Flat Connection	Ti 6Al 4V ELI
	FLAT6013	Titanium Gluing Ring Flat Connection	Ti 6Al 4V ELI
	FLAT6017	Long Laboratory Screw Flat Connectio	Ti 6Al 4V ELI
	FLAT6018	Centering Ring Flat Connection	Ti 6Al 4V ELI
	FLAT5512	Transfer Flat Connection	Ti 6Al 4V ELI

Ball Attachment - Overdentures

	MSD0051	SD Ball Attachment 0.5mm	Ti 6Al 4V ELI
	MSD0052	SD Ball Attachment 1mm	Ti 6Al 4V ELI
	MSD0053	SD Ball Attachment 2mm	Ti 6Al 4V ELI
	MSD0054	SD Ball Attachment 3mm	Ti 6Al 4V ELI
	MSD0055	SD Ball Attachment 4mm	Ti 6Al 4V ELI
	MSD0056	SD Ball Attachment 5mm	Ti 6Al 4V ELI
	NC01	Plastic Ball Cap White	Nylon
	NCP	Plastic Ball Cap Pink Soft	Nylon
	HOT	Titanium Ball Cap	Ti 6Al 4V ELI

Analogs

	MSD0007	SD Implant Analog	Ti 6Al 4V ELI
	MIH5740	Ball Attachment Abutment Replica	Stainless Steel

Transfers / Impression Coping

	MSD5109	SD Fix Konisch™ Slim Close Tray Transfer 3.6mmDx9mmL	Ti 6Al 4V ELI
	MSD5113	SD Fix Konisch™ Slim Close Tray Transfer 3.6mmDx13mmL	Ti 6Al 4V ELI
	MSD5209	SD Fix Konisch™ Close Tray Transfer 4.2mmDx9mmL	Ti 6Al 4V ELI
	MSD5213	SD Fix Konisch™ Close Tray Transfer 4.2mmDx13mmL	Ti 6Al 4V ELI
	MSD5610	SD Fix Konisch™ Slim Open Tray Transfer 3.6mmDx10mmL	Ti 6Al 4V ELI
	MSD5614	SD Fix Konisch™ Slim Open Tray Transfer 3.6mmDx14mmL	Ti 6Al 4V ELI
	MSD5510	SD Fix Konisch™ Open Tray Transfer 4.2mmDx10mmL	Ti 6Al 4V ELI
	MSD5514	SD Fix Konisch™ Open Tray Transfer 4.2mmDx14mmL	Ti 6Al 4V ELI

* Transfers can be tightened by Star Driver.

* Transfers will be chosen by length and width according to the planned prosthetic treatment and according to the interocclusal distance.

SD Retrieval Screw

	MSD3401	SD Retrieval Screw for Abutments	Ti 6Al 4V ELI
---	---------	----------------------------------	---------------

SD Screws

	MSD0073	SD Conical Implant Cover Screw	Ti 6Al 4V ELI
	MSD0001	SD Abutment Screw	Ti 6Al 4V ELI
	MSD0039	SD Retaining Screw for MU 1mm	Ti 6Al 4V ELI
	MSD0040	SD Retaining Screw for MU 2mm	Ti 6Al 4V ELI
	MSD0041	SD Retaining Screw for MU 3mm	Ti 6Al 4V ELI
	MSD3409	SD Retaining Screw for MU 4mm	Ti 6Al 4V ELI
	MSD3412	SD Retaining Screw for MU 5mm	Ti 6Al 4V ELI
	MSD3403	SD Retaining Screw for Angled MU	Ti 6Al 4V ELI
	MSD5010	SD Fix Konisch™ Open Tray Transfer Screw Short	Ti 6Al 4V ELI
	MSD5014	SD Fix Konisch™ Open Tray Transfer Screw Long	Ti 6Al 4V ELI
	MSD5126	SD Fix Konisch™ Close Tray Transfer Screw Short	Ti 6Al 4V ELI
	MSD5127	SD Fix Konisch™ Close Tray Transfer Screw Long	Ti 6Al 4V ELI

*All systems are provided with their compatible screws. However, they may also be ordered individually

Hand Drivers

	MSD0064	Star Hand Driver - Short	Stainless Steel
	MSD0004	Star Hand Driver - Long	Stainless Steel
	TMA-SHHS	TMA/FC Female Hand Hex Driver - Short	Stainless Steel
	TMA-SHHL	TMA/FC Female Hand Hex Driver - Long	Stainless Steel

Torque Drivers

	MSD0067	SD 2.2mm Hex Torque Driver-Short	Stainless Steel
	MSD0066	SD 2.2mm Hex Torque Driver-Long	Stainless Steel
	MSD0065	Star Torque Driver - Short	Stainless Steel
	MSD0005	Star Torque Driver - Long	Stainless Steel
	TDBA	Torque Driver for Ball Attachment	Stainless Steel
	TMA-SHT	TMA/FC Hex Torque Driver - Short	Stainless Steel
	TMA-SHTL	TMA/FC Hex Torque Driver - Long	Stainless Steel
	WRS	Wrench Ratchet Surgical Instruments Handle * use with MSD0068/ MSD0071/ MSD0074	Stainless Steel

Handpiece Drivers and Adaptors

	MSD0068	SD 2.25mm Hex Handpiece Adapter - Short	Stainless Steel
	MSD0071	SD 2.25mm Hex Handpiece Adapter - Medium	Stainless Steel
	MSD0074	SD 2.25mm Hex Handpiece Adapter - Long	Stainless Steel
	MSD0018	Star Handpiece Driver - Short	Stainless Steel
	MSD0019	Star Handpiece Driver - Long	Stainless Steel
	TMA-HAS	TMA/FC Hex Hand Piece Adapter - Short	Stainless Steel
	TMA-HAL	TMA/FC Hex Hand Piece Adapter - Long	Stainless Steel

Hand Drivers

	MSD0064	Star Hand Driver-Short	Stainless Steel
	MSD0004	Star Hand Driver - Longw	Stainless Steel
	TMA-SHHS	TMA/FC Female Hand Hex Driver - Short	Stainless Steel
	TMA-SHHL	TMA/FC Female Hand Hex Driver - Long	Stainless Steel

Torque Drivers

	MWD0067	WD 2.6mm Torque Driver-Short	Stainless Steel
	MWD0066	WD 2.6mm Torque Driver-Long	Stainless Steel
	MSD0065	Star Torque Driver-Short	Stainless Steel
	MSD0005	Star Torque Driver - Long	Stainless Steel
	TDBA	Torque Driver for Ball Attachment	Stainless Steel
	TMA-SHT	TMA/FC Hex Torque Driver - Short	Stainless Steel
	TMA-SHTL	TMA/FC Hex Torque Driver - Long	Stainless Steel
	WRS	Wrench Ratchet Surgical Instruments Handle	Stainless Steel

Handpiece Drivers and Adaptors

	MWD0068	WD 2.65mm Hex Handpiece Adapter - Short	Stainless Steel
	MWD0071	WD 2.65mm Hex Handpiece Adapter - Medium	Stainless Steel
	MWD0074	WD 2.65mm Hex Handpiece Adapter - Long	Stainless Steel
	MSD0018	Star Handpiece Driver - Short	Stainless Steel
	MSD0019	Star Handpiece Driver - Long	Stainless Steel
	TMA-HAS	TMA/FC Hex Hand Piece Adapter - Short	Stainless Steel
	TMA-HAL	TMA/FC Hex Hand Piece Adapter - Long	Stainless Steel

Fix-1™

One Piece Implant System

The Fix-1™ is a tapered core spiral implant.
It is as easy to use as regular crown and bridge dentistry.

Fix-1™ features:

- Minimally invasive surgery in narrow diameter implants
- Immediate function/loading
- Designed for narrow places
- Offers fast and easy flapless procedure
- Fix-1™ implants with narrow diameters are highly suitable replacements for laterals in the maxilla and incisors in the mandible

	L-10mm	L-11.5mm	L-13mm	L-15mm
Ø3.0mm	FO-3010	FO-3011	FO-3013	FO-3015
Ø3.3mm	FO-3310	FO-3311	FO-3313	FO-3315
Ø3.6mm	FO-3610	FO-3611	FO-3613	FO-3615
Ø4.2mm	FO-4210	FO-4211	FO-4213	FO-4215
Ø5.0mm	FO-5010	FO-5011	FO-5013	FO-5015

Packed with:



CA6025

One Piece Driver


Stainless Steel



AB/AE Surface
Treatment


Fix-1™ Ø 3.0

- ▶ AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

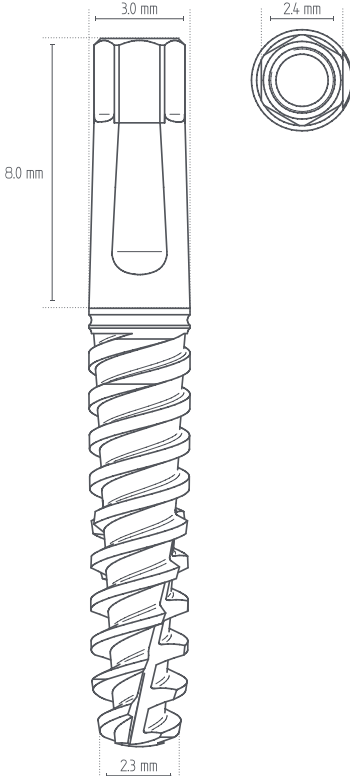Catalog No.	Length mm		
FO-3010	10		
FO-3011	11.5		
FO-3013	13		
FO-3015	15		

- ▶ Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. 2.0	1. 2.0 2. (2.8)	1. 2.0 2. 2.8

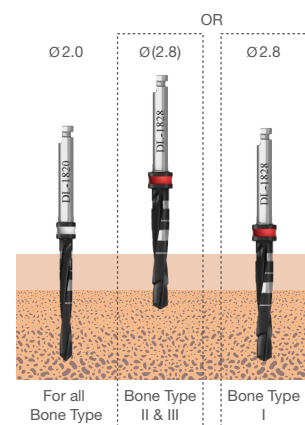


Note: (x.xx) drill to the depth of the cortex only. All measurements are in mm.

⚠ CAUTION: The drill preparation is up to 1mm longer than the implant.


Fix-1™ Ø 3.3

- ▶ AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

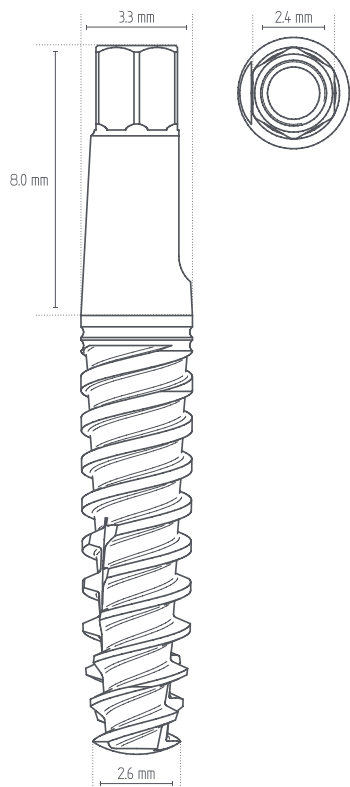Catalog No.	Length mm		
FO-3310	10		
FO-3311	11.5		
FO-3313	13		
FO-3315	15		

- ▶ Drill Sequence


BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. 2.0	1. 2.0 2. (2.8)	1. 2.0 2. 2.8

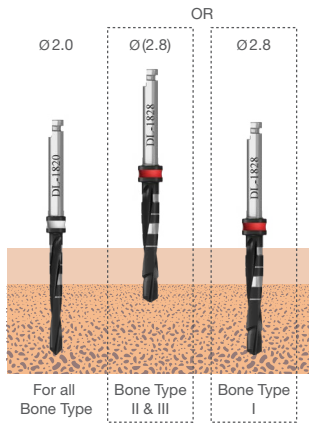

Note: (x.xx) drill to the depth of the cortex only. All measurements are in mm.
⚠ CAUTION: The drill preparation is up to 1mm longer than the implant.


Fix-1™ Ø 3.6

- ▶ AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

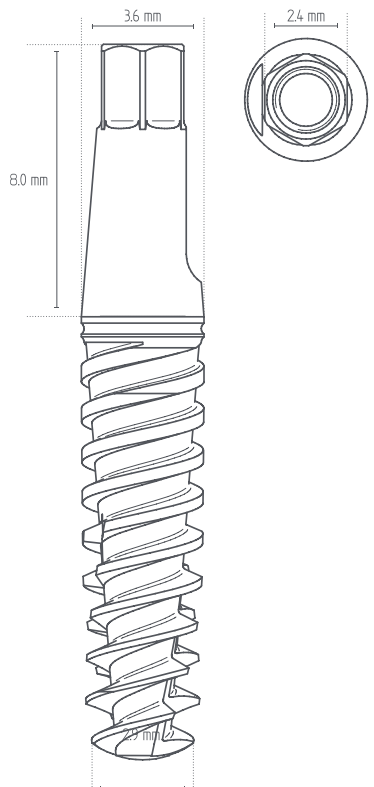Catalog No.	Length mm		
FO-3610	10		
FO-3611	11.5		
FO-3613	13		
FO-3615	15		

- ▶ Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step*	1. Tri-Step* 2. (3.2)	1. Tri-Step* 2. (3.2)

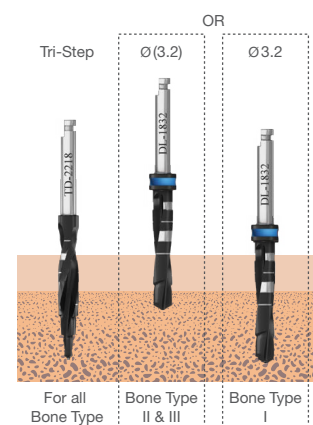

Note: (x.xx) drill to the depth of the cortex only. All measurements are in mm.




CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


Fix-1™ Ø 4.2

- ▶ AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

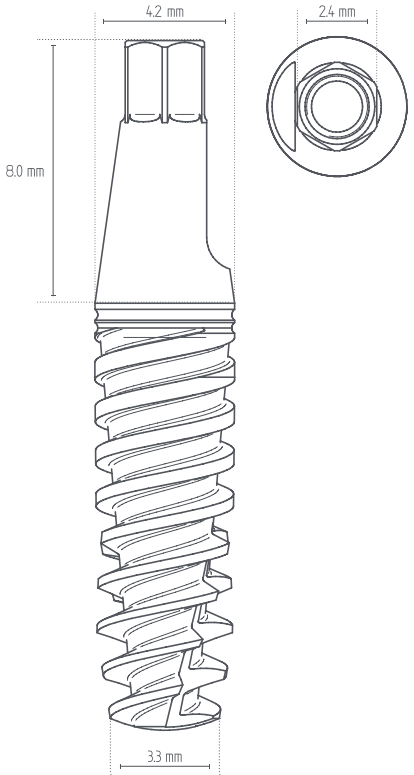Catalog No.	Length mm		
FO-4210	10		
FO-4211	11.5		
FO-4213	13		
FO-4215	15		

▶ Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step* 2. (3.6)	1. Tri-Step* 2. (3.6)	1. Tri-Step* 2. 3.6

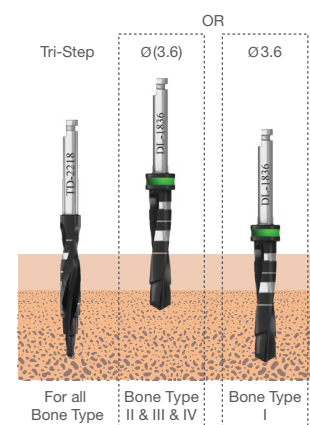

Note: (x.xx) drill to the depth of the cortex only. All measurements are in mm.



CAUTION: The drill preparation is up to 1mm longer than the implant.

* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.


Fix-1™ Ø 5.0

- ▶ AB/AE Surface
Titanium Alloy - Ti 6Al 4V ELI

Catalog No.	Length mm	
FO-5010	10	
FO-5011	11.5	
FO-5013	13	
FO-5015	15	

▶ Drill Sequence

BONE TYPE D-IV	BONE TYPE D-II & III	BONE TYPE D-I
1. Tri-Step* 2. 3.6 3. (4.2)	1. Tri-Step* 2. 3.6 3. 4.2 4. (4.6)	1. Tri-Step* 2. 3.6 3. 4.2 4. 4.6

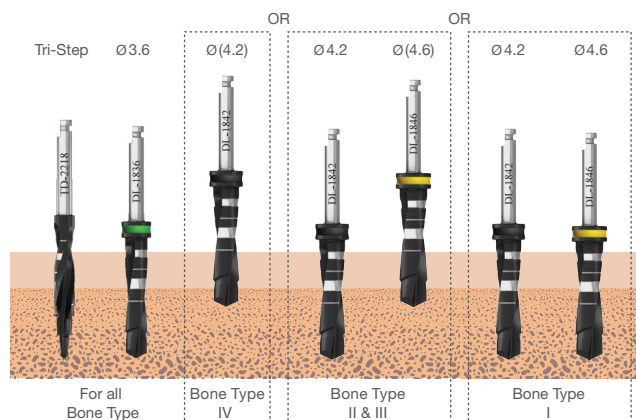

Note: (x.xx) drill to the depth of the cortex only.
All measurements are in mm.

⚠ CAUTION: The drill preparation is up to 1mm longer than the implant.
* For initial drilling, you may use 2, 2.8 and 3.2 drills in sequence instead of the Tri-Step drill.

CAD CAM






CAD CAM IH

Regular Standard
Internal Hex

Ø 3.3 Ø 3.5 Ø 3.75
Ø 4.2 Ø 5 Ø 6


FixTite™ - S | Mars™ | Fix-a-Dent™

Scan body | For Intra-oral and Desktop Scanner


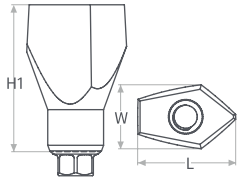






IH Scan Body 7mm	IH Scan Body 10mm
	
Catalog No. MIH0088	Catalog No. MIH0089
Retaining Screw Included SRAS	Retaining Screw Included MIH3425
Material TI 6Al 4V ELI	Material TI 6Al 4V ELI
Dimensions H1: 7 mm W: 4.3 mm L: 6.7 mm	Dimensions H1: 10 mm W: 4.4 mm L: 6.8 mm

Ti Blanks | For Medentika PreFace® Holder

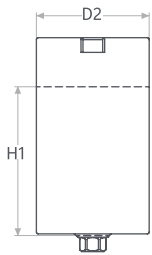


IH Pre Milled MDTK Ti Blank Ø11.5mm	IH Pre Milled MDTK Ti Blank Ø15.8mm
	
Catalog No. MIH1411	Catalog No. MIH1415
Retaining Screw Included SRAS2	Retaining Screw Included SRAS2
Material TI 6Al 4V ELI	Material TI 6Al 4V ELI
Dimensions H1: 15.2mm D2: Ø11.5 mm	Dimensions H1: 15.2 mm D2: Ø15.8 mm
Intended Use customized abutment for single unit	Intended Use customized abutment for single unit

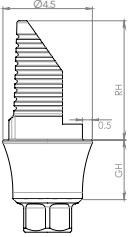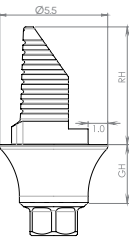Ti Bases		RH	GH	Catalog No.
	 IH T-Base, Esthetic	4.0	1.0	MIH5142
		4.0	2.0	MIH5148
		4.0	3.0	MIH5154
		6.0	1.0	MIH5166
		6.0	2.0	MIH5172
		6.0	3.0	MIH5178
		8.0	1.0	MIH5190
		8.0	2.0	MIH5196
		8.0	3.0	MIH5202
	 IH T-Base, Esthetic Wide	4.0	1.0	MIH5143
		4.0	2.0	MIH5149
		4.0	3.0	MIH5155
		6.0	1.0	MIH5167
		6.0	2.0	MIH5173
		6.0	3.0	MIH5179
		8.0	1.0	MIH5191
		8.0	2.0	MIH5197
		8.0	3.0	MIH5203

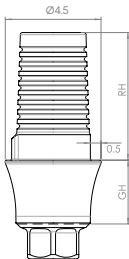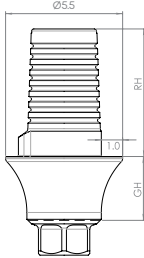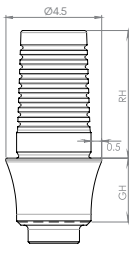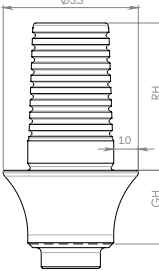Ti Bases	Description	RH	GH	Catalog No.
	 IH T-Base, Anti rotation	4.0	1.0	MIH5138
		4.0	2.0	MIH5144
		4.0	3.0	MIH5150
		4.0	4.0	MIH5156
		6.0	1.0	MIH5162
		6.0	2.0	MIH5168
		6.0	3.0	MIH5174
		6.0	4.0	MIH5180
		8.0	1.0	MIH5186
		8.0	2.0	MIH5192
		8.0	3.0	MIH5198
		8.0	4.0	MIH5204
	 IH T-Base, Anti rotation Wide	4.0	1.0	MIH5139
		4.0	2.0	MIH5145
		4.0	3.0	MIH5151
		4.0	4.0	MIH5157
		6.0	1.0	MIH5163
		6.0	2.0	MIH5169
		6.0	3.0	MIH5175
		6.0	4.0	MIH5181
		8.0	1.0	MIH5187
		8.0	2.0	MIH5193
		8.0	3.0	MIH5199
		8.0	4.0	MIH5205
	 IH T-Base, Free rotation	4.0	1.0	MIH5140
		4.0	2.0	MIH5146
		4.0	3.0	MIH5152
		4.0	4.0	MIH5158
		6.0	1.0	MIH5164
		6.0	2.0	MIH5170
		6.0	3.0	MIH5176
		6.0	4.0	MIH5182
		8.0	1.0	MIH5188
		8.0	2.0	MIH5194
		8.0	3.0	MIH5200
		8.0	4.0	MIH5206
	 IH T-Base, Free rotation Wide	4.0	1.0	MIH5141
		4.0	2.0	MIH5147
		4.0	3.0	MIH5153
		4.0	4.0	MIH5159
		6.0	1.0	MIH5165
		6.0	2.0	MIH5171
		6.0	3.0	MIH5177
		6.0	4.0	MIH5183
		8.0	1.0	MIH5189
		8.0	2.0	MIH5195
		8.0	3.0	MIH5201
		8.0	4.0	MIH5207

Scan body | For Intra-oral and Desktop Scanner

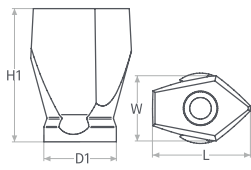

Catalog No.
Retaining Screw Included
Material
Dimensions


MU Scan Body 6mm

MMU0006
MMU0001
TI 6Al 4V ELI
H1: 6 mm D1: Ø4.9 mm W: 4.3 mm L: 6.7 mm

MU Scan Body 9mm

MMU0009
MMU0001
TI 6Al 4V ELI
H1: 9 mm D1: Ø4.9 mm W: 4.5 mm L: 6.8 mm


Cementing Cone & Cylinder

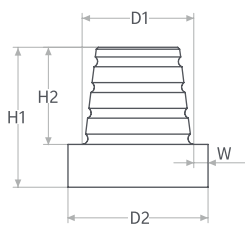

Catalog No.
Retaining Screw Included
Material
Dimensions

MU Cementing Cone


MMU0010
MMU0001
TI 6Al 4V ELI
H1: 4.9 mm H2: 3.4 mm D1: Ø3.8 mm D2: Ø4.9 mm W: 0.5 mm

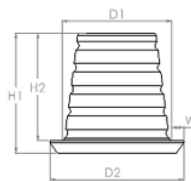

Catalog No.
Retaining Screw Included
Material
Dimensions

MU Cementing Con Low Profile

MMU0011
MMU0001
TI 6Al 4V ELI
H1: 4.9 mm H2: 4.45 mm D1: Ø4.5 mm D2: Ø5.5 mm W: 0.5 mm

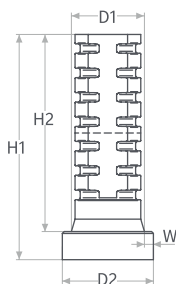

Catalog No.
Retaining Screw Included
Material
Dimensions

Titanium Temporary Copying

TMA-TTC
TPS
TI 6Al 4V ELI
H1: 12 mm H2: 10.5 mm D1: Ø3.9 mm D2: Ø4.9 mm W: 0.5 mm




CAD CAM

ND Fix Konisch™
Conical-Hex Connection


Ø 3

Scan body | For Intra-oral and Desktop Scanner



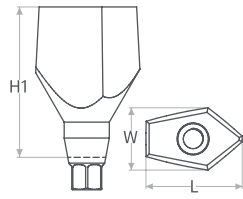


	ND Fix Konisch™ Scan Body 7mm	ND Fix Konisch™ Scan Body 10mm
Catalog No.	MND0088	MND0089
Retaining Screw Included	MND0001	MND0001
Material	TI 6Al 4V ELI	TI 6Al 4V ELI
Dimensions	H1: 7 mm W: 4.2 mm L: 6.7 mm	H1: 10 mm W: 4.5 mm L: 6.9 mm

Ti Blanks | For Medentika PreFace® Holder

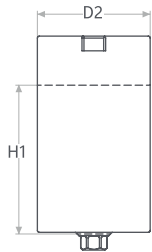

	ND Fix Konisch™ Pre Milled MDTK Ti Blank Ø 11.5mm		ND Fix Konisch™ Pre Milled MDTK Ti Blank Ø 11.5mm
Catalog No.	MND1111	Catalog No.	MND1115
Retaining Screw Included	MND0001	Retaining Screw Included	MND0001
Material	TI 6Al 4V ELI	Material	TI 6Al 4V ELI
Dimensions	H1: 15.4 mm D2: Ø 11.5 mm	Dimensions	H1: 15.4 mm D2: Ø 15.8 mm
Intended Use	customized abutment for single unit	Intended Use	customized abutment for single unit

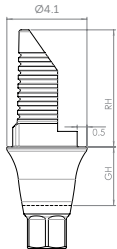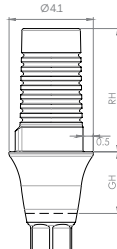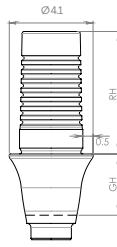Ti Bases		RH	GH	Catalog No.
	 ND T-base, Esthetic	4.0	1.5	MND5002
		4.0	2.0	MND5005
		4.0	3.0	MND5008
		6.0	1.5	MND5013
		6.0	2.0	MND5017
		6.0	3.0	MND5020
		8.0	1.5	MND5026
		8.0	2.0	MND5029
	 ND T-base, Anti rotation	4.0	1.5	MND5000
		4.0	2.0	MND5003
		4.0	3.0	MND5006
		4.0	4.0	MND5009
		6.0	1.5	MND5011
		6.0	2.0	MND5015
		6.0	3.0	MND5018
		6.0	4.0	MND5021
		8.0	1.5	MND5024
		8.0	2.0	MND5027
	 ND T-base, Free rotation	4.0	1.5	MND5001
		4.0	2.0	MND5004
		4.0	3.0	MND5007
		4.0	4.0	MND5010
		6.0	1.5	MND5012
		6.0	2.0	MND5016
		6.0	3.0	MND5019
		6.0	4.0	MND5022
		8.0	1.5	MND5025
		8.0	2.0	MND5028




CAD CAM

Fix Konisch™ PRO

Standard/Wide
Ø3.5 Ø4.3 Ø5.0

Standard Diameter for size ridges and spaces

Scan body | For Intra-oral and Desktop Scanner

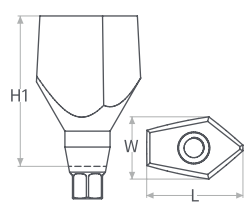

SD Fix Konisch™
Scan Body 7mm



SD Fix Konisch™
Scan Body 10mm



Catalog No.
Retaining Screw Included
Material
Dimensions

MSD0088

MSD0001

TI 6Al 4V ELI

H1: 7 mm
W: 4.3 mm
L: 6.7 mm


MSD0089

MSD0001


TI 6Al 4V ELI

H1: 10 mm
W: 4.5 mm
L: 6.9 mm

Ti Blanks | For Medentika PreFace® Holder

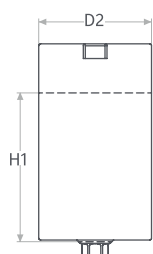

SD Fix Konisch™
Pre Milled
MDTK Ti Blank
Ø11.5mm

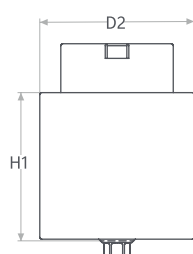

SD Fix Konisch™
Pre Milled
MDTK Ti Blank
Ø15.8mm



Catalog No.
Retaining Screw Included
Material
Dimensions
Intended Use

MSD1211

MSD0001

TI 6Al 4V ELI

H1: 15.2 mm
D2: Ø 11.5 mm

customized abutment
for single unit

Catalog No.
Retaining Screw Included
Material
Dimensions
Intended Use




MSD1215

MSD0001

TI 6Al 4V ELI

H1: 15.2 mm
D2: Ø 15.8 mm

customized abutment
for single unit

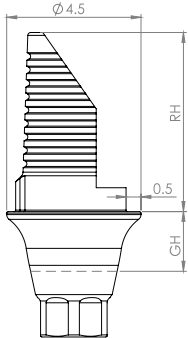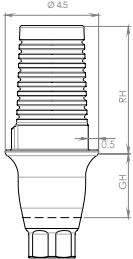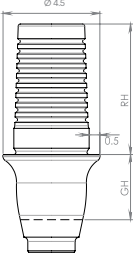Ti Bases		RH	GH	Catalog No.
	 PRO T-base, Esthetic	4.0	1.5	MSD5032
		4.0	2.0	MSD5035
		4.0	3.0	MSD5038
		6.0	1.5	MSD5044
		6.0	2.0	MSD5047
		6.0	3.0	MSD5050
		8.0	1.5	MSD5056
		8.0	2.0	MSD5059
		8.0	3.0	MSD5062
	 PRO T-base, Anti rotation	4.0	1.0	MSD5030
		4.0	2.0	MSD5033
		4.0	3.0	MSD5036
		4.0	4.0	MSD5039
		6.0	1.0	MSD5042
		6.0	2.0	MSD5045
		6.0	3.0	MSD5048
		6.0	4.0	MSD5051
		8.0	1.0	MSD5054
		8.0	2.0	MSD5057
		8.0	3.0	MSD5060
		8.0	4.0	MSD5063
	 PRO T-base, Free rotation	4.0	1.0	MSD5031
		4.0	2.0	MSD5034
		4.0	3.0	MSD5037
		4.0	4.0	MSD5040
		6.0	1.0	MSD5043
		6.0	2.0	MSD5046
		6.0	3.0	MSD5049
		6.0	4.0	MSD5052
		8.0	1.0	MSD5055
		8.0	2.0	MSD5058
		8.0	3.0	MSD5061
		8.0	4.0	MSD5064



Accessories

Screws, Tools, Analogs & Replica

Screws

Fix Konisch™

	ND Abutment Screw	SD Abutment Screw
		
Catalog No.	MND0001	MSD0001
Material	Ti 6Al 4V ELI	Ti 6Al 4V ELI
Intended use	For ND Scan Body, Bases and Ti Blanks	For SD Scan Body, Bases and Ti Blanks

IH Regular Standard Internal Hex

	IH Prosthetic Screw	IH Screw For Closed Tray Transfer	IH Short Prosthetic Screw
			
Catalog No.	SRAS	MIH3425	SRAS2
Material	Ti 6Al 4V ELI	Ti 6Al 4V ELI	Ti 6Al 4V ELI
Intended Use	For IH Scan Body 7mm	For IH Scan Body 10mm	For IH Ti Bases and Ti Blanks

MU - Multi Unit

	MU Scan Body Screw	MU Prosthetic Screw
		
Catalog No.	MMU0001	TPS
Material	Ti 6Al 4V ELI	Ti 6Al 4V ELI
Intended Use	For MU Scan Body and Cementing Cone	For MU Cylinder and MU direct option

Analogs & Replica*

Digital analog for Fix Konisch™

	ND Digital Model Analog	SD Digital Model Analog
		
Catalog No.	MND1000	MSD1000
Material	TI 6Al 4V ELI	TI 6Al 4V ELI

Digital Analog for IH and MU

	IH Digital Model Analog	MU Digital Model abutment Replica
		
Catalog No.	MIH1000	MIH1001
Material	TI 6Al 4V ELI	TI 6Al 4V ELI

Tools

	Star Hand Driver- Short	1.27 Hex Hand Driver- Short
		
Catalog No.	MSD0064	SHHS
	Star Hand Driver- Long	1.27 Hex Hand Driver- Medium
Catalog No.	MSD0004	SHHM
Intended Use	Compatible with MU and Fix Konisch™ Platforms: ND, SD, WD	Compatible with IH System: FixTite™-S, Mars™, Fix-a-Dent™

*You may also use any of the following:

SHHL	MSD0065
SHT	MSD0005
SHTL	MSD0018
MIH9030	MSD0019
MIH9035	FIA

Kits & Tools

Standard Internal Hex Surgical Kit Premium



Also available with short drills
SURPSS Standard Internal Hex Plastic Surgical Kit - Starter Short Drills

SURPP: Standard Internal Hex Plastic Surgical Kit - Premium**DRILLS:**

MTD1	Tip Drill
2218TDM	Three Step Drill - Long
1820DM	Pilot Drill 2mmD x 8-18mmL
1828DM	Twist Step Drill 2.3/2.8mmD x 8-18mmL
1832DM	Twist Step Drill 2.8/3.2mmD x 8-18mmL
1839DM	Twist Step Drill 3.2/3.9mmD x 8-18mmL
1836DM	Twist Step Drill 3.2/3.6mmD x 8-18mmL
1842DM	Twist Step Drill 3.8/4.2mmD x 8-18mmL
1846DM	Twist Step Drill 4.2/4.6mmD x 8-18mmL
1852DM	Twist Step Drill 4.8/5.2mmD x 8-18mmL
1511TDM	Three Step Drill - Short
1320DM	Pilot Drill 2mmD x 6-13mmL
1328DM	Twist Step Drill 2.3/2.8mmDx6-13mmL
1332DM	Twist Step Drill 2.8/3.2mmDx6-13mmL
1336DM	Twist Step Drill 3.2/3.6mmDx6-13mmL
1339DM	Twist Step Drill 3.2/3.9mmD x 6-13mmL
1342DM	Twist Step Drill 3.8/4.2mmDx6-13mmL
1346DM	Twist Step Drill 4.2/4.6mmDx6-13mmL
1352DM	Twist Step Drill 4.8/5.2mmDx6-13mmL
1356DM	Twist Step Drill 5.2/5.6mmDx6-13mmL

PARALLEL PINS:

PAPS	Parallel Pin 11 (qty. 4)
PAPL	Parallel Pin 16 (qty. 2)

TOOLS:

SHHA9021	Hex Handpiece Adapter - Short
MHHA9026	Hex Handpiece Adapter - Medium
LHHA9026	Hex Handpiece Adapter - Long

HEX HAND DRIVERS:

SHHS	1.27 Hex Hand Driver - Short
SHHM	1.27 Hex Hand Driver - Medium

HEX TORQUE DRIVERS:

SHT	1.27 Hex Torque Driver - Short
SHTL	1.27 Hex Torque Driver - Long

DRILL EXTENSION:

ITHS	Drill Extension
SRS	Surgical Torque Ratchet 35-100 Ncm
WRS	Wrench Ratchet Surgical Instruments Handle
DP6180	Depth Probe

SURPA: Standard Internal Hex Plastic Surgical Kit - Advanced**DRILLS:**

1820DM	Pilot Drill 2mmD x 8-18mmL
1828DM	Twist Step Drill 2.3/2.8mmD x 8-18mmL
1832DM	Twist Step Drill 2.8/3.2mmD x 8-18mmL
1836DM	Twist Step Drill 3.2/3.6mmD x 8-18mmL
1842DM	Twist Step Drill 3.8/4.2mmD x 8-18mmL
1846DM	Twist Step Drill 4.2/4.6mmD x 8-18mmL

PARALLEL PINS:

PAPS	Parallel Pin 11 (qty. 4)
PAPL	Parallel Pin 16 (qty. 2)
WRS	Wrench Ratchet Surgical Instruments Handle

TOOLS:

SHHA9021	RS Hex Handpiece Adapter Short
MHHA9026	RS Hex Handpiece Adapter Medium
LHHA9029	RS Hex Hand Piece Adapter Long

HEX HAND DRIVERS:

SHHS	1.27 Hex Hand Driver - Short
SHHM	1.27 Hex Hand Driver - Medium

HEX TORQUE DRIVERS:

SHT	1.27 Hex Torque Driver - Short
SHTL	1.27 Hex Torque Driver - Long

DRILL EXTENSION:

ITHS	Drill Extension
SRS	Surgical Torque Ratchet 35-100 Ncm
DP6180	Depth Probe

Also available with short drills

SURPAS - Standard Internal Hex Plastic Surgical Kit - Advanced Short Drills

SURPS: Standard Internal Hex Plastic Surgical Kit - starter**DRILLS:**

1820DM	Pilot Drill 2mmD x 8-18mmL
1828DM	Twist Step Drill 2.3/2.8mmD x 8-18mmL
1832DM	Twist Step Drill 2.8/3.2mmD x 8-18mmL
1836DM	Twist Step Drill 3.2/3.6mmD x 8-18mmL
1842DM	Twist Step Drill 3.8/4.2mmD x 8-18mmL
1846DM	Twist Step Drill 4.2/4.6mmD x 8-18mmL

TOOLS:

SHHA9021	RS Hex Handpiece Adapter Short
MHHA9026	RS Hex Handpiece Adapter Medium
LHHA9029	RS Hex Hand Piece Adapter Long

HEX HAND DRIVERS:

SHHS	1.27 Hex Hand Driver - Short
SHHM	1.27 Hex Hand Driver - Medium

SR	Wrench Ratchet
WRS	Wrench Ratchet Surgical Instruments Handle

Fix Konisch™ Surgical Kit Premium



FK6300: Fix Konisch™ Surgical Kit - Premium

SRS	Surgical Torque Ratchet 35-100 Ncm	1342DM	Twist Step Drill 3.8/4.2mmD x 6-13mmL
DP6180	Depth Probe	1346DM	Twist Step Drill 4.2/4.6mmD x 6-13mmL
PAPS	Parallel Pin 11	ITHS	Drill Extension
PAPL	Parallel Pin 16	MSD0064	Star Hand Driver-Short
MTD1	Tip Drill	MSD0004	Star Hand Driver - Long
1511TDM	Three Step Drill-Short	MSD0065	Star Torque Driver - Short
2218TDM	Three Step Drill-Long	MSD0005	Star Torque Driver - Long
1820DM	Pilot Drill 2mmD x 8-18mmL	MND0068	ND 2mm Hex Handpiece Adapter - Short
1828DM	Twist Step Drill 2.3/2.8mmD x 8-18mmL	MND0071	ND 2mm Hex Handpiece Adapter - Medium
1832DM	Twist Step Drill 2.8/3.2mmD x 8-18mmL	MND0074	ND 2mm Hex Handpiece Adapter - Long
1836DM	Twist Step Drill 3.2/3.6mmD x 8-18mmL	MSD0068	SD 2.25mm Hex Handpiece Adapter - Short
1842DM	Twist Step Drill 3.8/4.2mmD x 8-18mmL	MSD0071	SD 2.25mm Hex Handpiece Adapter - Medium
1846DM	Twist Step Drill 4.2/4.6mmD x 8-18mmL	WRS	Wrench Ratchet Surgical Instruments Handl
1320DM	Pilot Drill 2mmD x 6-13mmL		
1328DM	Twist Step Drill 2.3/2.8mmD x 6-13mmL		
1332DM	Twist Step Drill 2.8/3.2mmD x 6-13mmL		
1336DM	Twist Step Drill 3.2/3.6mmD x 6-13mmL		

FK6302: Fix Konisch™ Surgical Kit - Advanced

SRS	Surgical Torque Ratchet 35-100 Ncm	MSD0004	Star Hand Driver - Long
DP6180	Depth Probe	MSD0065	Star Torque Driver - Short
PAPS	Parallel Pin 11	MSD0005	Star Torque Driver - Long
PAPL	Parallel Pin 16	MND0068	ND 2mm Hex Handpiece Adapter - Short
1820DM	Pilot Drill 2mmD x 8-18mmL	MND0071	ND 2mm Hex Handpiece Adapter - Medium
1828DM	Twist Step Drill 2.3/2.8mmD x 8-18mmL	MND0074	ND 2mm Hex Handpiece Adapter - Long
1832DM	Twist Step Drill 2.8/3.2mmD x 8-18mmL	MSD0068	SD 2.25mm Hex Handpiece Adapter - Short
1836DM	Twist Step Drill 3.2/3.6mmD x 8-18mmL	MSD0071	SD 2.25mm Hex Handpiece Adapter - Medium
1842DM	Twist Step Drill 3.8/4.2mmD x 8-18mmL	MSD0074	SD 2.25mm Hex Handpiece Adapter - Long
1846DM	Twist Step Drill 4.2/4.6mmD x 8-18mm	WRS	Depth Wrench Ratchet Surgical Instruments Handle
ITHS	Drill Extension		
MSD0064	Star Hand Driver-Short		

Tools

Tools

SRS



Surgical Torque Ratchet 35-100 Ncm

Stainless Steel

PTR



Prosthetic Torque Ratchet 15-35 Ncm

Stainless Steel

SR



Wrench Ratchet

Stainless Steel

SDH1



Surgical Driver

Stainless Steel

DP6180



Depth Probe

Stainless Steel

PAPS



Parallel Pin 11

Stainless Steel

PAPL



Parallel Pin 16

Stainless Steel

Osteotomes

SOS-H



Osteotome Handle

Aluminium + Stainless Steel

SOS20



Straight Convex Osteotome 2.0mmD

Ti 6Al 4V ELI

SOS325



Straight Convex Osteotome 3.25mmD

Ti 6Al 4V ELI

SOS28



Straight Convex Osteotome 2.8mmD

Ti 6Al 4V ELI

SOS365



Straight Convex Osteotome 3.65mmD

Ti 6Al 4V ELI

Standard Internal Hex Surgical Tools

	SMTS	2.4 Hex Torque Driver - Short	Stainless
	SMTL	2.4 Hex Torque Driver - Long	Stainless
	LHHA9029	RS Hex Hand Piece Adapter Long	Stainless Steel
	MHHA9026	RS Hex Handpiece Adapter Medium	Stainless Steel
	SHHA9021	RS Hex Handpiece Adapter Short	Stainless Steel

Standard Internal Hex Prosthetic Tools

	SHT	1.27 Hex Torque Driver - Short	Stainless Steel
	SHTL	1.27 Hex Torque Driver - Long	Stainless Steel
	SHHM	1.27 Hex Hand Driver - Medium	Stainless Steel
	SHHS	1.27 Hex Hand Driver - Short	Stainless Steel
	SHHL	1.27 Hex Hand Driver - Long	Stainless Steel
	TDBA	Torque Driver for Ball Attacment	Stainless Steel
	TMA-HAS	TMA /FC Hex Handpiece Adapter - Short	Stainless Steel
	TMA-HAL	TMA /FC Hex Handpiece Adapter - Long	Stainless Steel
	TMA-SHT	TMA /FC Hex Torque Driver - Short	Stainless Steel
	TMA-SHTL	TMA /FC Hex Torque Driver - Long	Stainless Steel
	TMA-SHHS	TMA / FC Female Hand Hex Driver - Short	Stainless Steel
	TMA-SHHL	TMA / FC Female Hand Hex Driver - Long	Stainless Steel

Surgical Drills

Two Step DLC Drills - Long

	1820DM	Pilot Drill 2mmD x 8-18mmL White	White	Stainless Steel & DLC
	1828DM	Twist Step Drill 2.3/2.8mmD x 8-18mmL	Red	Stainless Steel & DLC
	1832DM	Twist Step Drill 2.8/3.2mmD x 8-18mmL	Blue	Stainless Steel & DLC
	1836DM	Twist Step Drill 3.2/3.6mmD x 8-18mmL	Green	Stainless Steel & DLC
	1839DM	Twist Step Drill 3.2/3.9mmD x 8-18mmL	Pink	Stainless Steel & DLC
	1842DM	Twist Step Drill 3.8/4.2mmD x 8-18mmL	Black	Stainless Steel & DLC
	1846DM	Twist Step Drill 4.2/4.6mmD x 8-18mmL	Yellow	Stainless Steel & DLC
	1852DM	Twist Step Drill 4.8/5.2mmD x 8-18mmL	Brown	Stainless Steel & DLC
	1856DM	Twist Step Drill 5.2/5.6mmD x 8-18mmL	Grey	Stainless Steel & DLC

Two Step DLC Drills - Short

	1320DM	Pilot Drill 2mmD x 6-13mmL	White	Stainless Steel & DLC
	1328DM	Twist Step Drill 2.3/2.8mmD x 6-13mmL	Red	Stainless Steel & DLC
	1332DM	Twist Step Drill 2.8/3.2mmD x 6-13mmL	Blue	Stainless Steel & DLC
	1336DM	Twist Step Drill 3.2/3.6mmD x 6-13mmL	Green	Stainless Steel & DLC
	1339DM	Twist Step Drill 3.2/3.9mmD x 6-13mmL	Pink	Stainless Steel & DLC
	1342DM	Twist Step Drill 3.8/4.2mmD x 6-13mmL	Black	Stainless Steel & DLC
	1346DM	Twist Step Drill 4.2/4.6mmD x 6-13mmL	Yellow	Stainless Steel & DLC
	1352DM	Twist Step Drill 4.8/5.2mmD x 6-13mmL	Brown	Stainless Steel & DLC
	1356DM	Twist Step Drill 5.2/5.6mmD x 6-13mmL	Grey	Stainless Steel & DLC

Three Step DLC Drills

	1511TDM	Three Step Drill-Short	Stainless Steel & DLC
	2218TDM	Three Step Drill-Long	Stainless Steel & DLC

Tip Drill

	MTD1	Tip Drill	Stainless Steel
---	------	-----------	-----------------

Drill Extension

	ITHS	Drill Extension	Stainless Steel
---	------	-----------------	-----------------

Drill Stoppers

Long Drill - Regular



Each Drill Stop is marked with width and length.
 Note: Some stops are suitable for long and short drill according to the following markings.
 L: Drill length for long drill (18 mm)
 S: Drill length for short drill (13 mm)

Short Drill 13mm

L	2.0mmD	2.8mmD	3.2mmD	3.6mmD	4.2mmD	4.6mmD	5.2mmD	5.6mmD
6mm	1320L6SM	1328L6SM	1332L6SM	1336L6SM	1342L6SM	1346L6SM	1352L6SM	1356L6SM
8mm	1820L13SM	1828L13SM	1832L13SM	1836L13SM	1842L13SM	1846L13SM	1852L13SM	1856L13SM
10mm	1820L15SM	1828L15SM	1832L15SM	1836L15SM	1842L15SM	1846L15SM	1852L15SM	1856L15SM
11.5mm	1320L11SM	1328L11SM	1332L11SM	1336L11SM	1342L11SM	1346L11SM	1352L11SM	1356L11SM

Long Drill 18mm

L	2.0mmD	2.8mmD	3.2mmD	3.6mmD	4.2mmD	4.6mmD	5.2mmD	5.6mmD
6mm	1820L6SM	1828L6SM	1832L6SM	1836L6SM	1842L6SM	1846L6SM	1852L6SM	1856L6SM
8mm	1820L8SM	1828L8SM	1832L8SM	1836L8SM	1842L8SM	1846L8SM	1852L8SM	1856L8SM
10mm	1820L10SM	1828L10SM	1832L10SM	1836L10SM	1842L10SM	1846L10SM	1852L10SM	1856L10SM
11.5mm	1820L11SM	1828L11SM	1832L11SM	1836L11SM	1842L11SM	1846L11SM	1852L11SM	1856L11SM
13mm	1820L13SM	1828L13SM	1832L13SM	1836L13SM	1842L13SM	1846L13SM	1852L13SM	1856L13SM
15mm	1820L15SM	1828L15SM	1832L15SM	1836L15SM	1842L15SM	1846L15SM	1852L15SM	1856L15SM
16mm	1820L16SM	1828L16SM	1832L16SM	1836L16SM	1842L16SM	1846L16SM	1852L16SM	1856L16SM

Drill Stop Kits	
S1001M	Long Drill Stop Kit - Regular
S1002M	Long Drill Stop Kit - Wide
S1003M	Short Drill Stop Kit - Regular
S1004M	Short Drill Stop Kit - Wide

Guided Surgery Kit



Tools

	Anchoring Drill S-2.0, 2.0/1.5mmDx19mmL	MGSD1915
	Anchoring Drill S-2.0, 2.0/1.7mmDx19mmL	MGSD1917
	Handpiece Guide Fixation Screw Driver	MGST0028
	Tissue Punch 5.0mmD, S-5.5	MGSD5055
	Drill Extension	ITHS
	Screw Fixture Mount	MGST0029
	Handpiece Fixture Mount	MGST0032
	Grip Fixture Mount	MGST0034
	Manual Socket L10.5	MGST0023
	Manual Socket L14.5	MGST0024
	Handpiece Socket L10.5	MGST0021
	Handpiece Socket L14.5	MGST0022
	Extraction Tool	MGST0035
	1.27 Hex Hand Driver - Medium	SHHM
	Guide Fixation Screw Driver	MGST0027
	Torque Ratchet 35-100 Ncm 4mm Hex	MGST0011

Drills		
	Pilot Drill S-5.5, 2.5/1.8mmDx6.0mmL	MGSD0625
	Drill S-5.5, 2.5/1.8mmDx16.0mmL	MGSD1625
	Drill S-5.5, 2.5/1.8mmDx13.0mmL	MGSD1325
	Drill S-5.5, 2.5/1.8mmDx11.5mmL	MGSD1125
	Drill S-5.5, 2.5/1.8mmDx10.0mmL	MGSD1025
	Drill S-5.5, 2.5/1.8mmDx8.0mmL	MGSD0825
	Drill S-5.5, 3.2/2.5mmDx16.0mmL	MGSD1632
	Drill S-5.5, 3.2/2.5mmDx13.0mmL	MGSD1332
	Drill S-5.5, 3.2/2.5mmDx11.5mmL	MGSD1132
	Drill S-5.5, 3.2/2.5mmDx10.0mmL	MGSD1032
	Drill S-5.5, 3.2/2.5mmDx8.0mmL	MGSD0832
	Drill S-5.5, 3.2/2.5mmDx6.25mmL	MGSD0632
	Drill S-5.5, 3.85/3.2mmDx16.0mmL	MGSD1638
	Drill S-5.5, 3.85/3.2mmDx13.0mmL	MGSD1338
	Drill S-5.5, 3.85/3.2mmDx11.5mmL	MGSD1138
	Drill S-5.5, 3.85/3.2mmDx10.0mmL	MGSD1038
	Drill S-5.5, 3.85/3.2mmDx8.0mmL	MGSD0838
	Drill S-5.5, 3.85/3.2mmDx6.25mmL	MGSD0638
	Drill S-5.5, 4.6/3.85mmDx16.0mmL	MGSD1646
	Drill S-5.5, 4.6/3.85mmDx13.0mmL	MGSD1346
	Drill S-5.5, 4.6/3.85mmDx11.5mmL	MGSD1146
	Drill S-5.5, 4.6/3.85mmDx10.0mmL	MGSD1046
	Drill S-5.5, 4.6/3.85mmDx8.0mmL	MGSD0846
	Drill S-5.5, 4.6/3.85mmDx6.25mmL	MGSD0646

Additional Products Sold separately

18mm Drills			
	Drill S-5.5, 2.5/1.8mmDx18.0mmL		MGSD1825
	Drill S-5.5, 3.2/2.5mmDx18.0mmL		MGSD1832
	Drill S-5.5, 3.85/3.2mmDx18.0mmL		MGSD1838
	Drill S-5.5, 4.6/3.85mmDx18.0mmL		MGSD1846

Accessories			
	Guided Vertical fixation screw S-5.5	x5 pack	MGST0007
	2.0D 19mmL Anchorage Screw	x5 pack	MGSD1920
	Fixation Screw 2.0mmD Sleeve	x5 pack	MGST0020
	Regular/ Wide Implant 5.5mmD Sleeve	x5 pack	MGST0055



www.medigma-implants.com

 +49 175 8726995

 Siemensstrasse 10 D-78564 Wehingen, Germany

 office@medigma-implants.com